

Electronic Signature Information

Name	DOC0896460
Revision	1
Type	Controlled Document
Title	Dicom Conformance Statement MammoWorkstation 4.6.0
Originator	100037834_robort__lopuszanski

Name	Reason For Change	File Size (Bytes)
Conformance Statement MammoWorkstation.doc	MammoWorkstation DICOM Conformance Statement	1413632

Route	Signer	Status	Comments	Completion Date
R-4418842	100037834_robort__lopuszanski	Approve	Approved as LSD	20 Dec 2010 11:11:40 GMT
R-4418842	100040664_clemens__hermann	Approve	approved as engineering manager	20 Dec 2010 12:16:40 GMT
R-4418944	100037834_robort__lopuszanski	Approve	Approved 2nd route as LSD	20 Dec 2010 12:40:52 GMT
R-4418944	100040664_clemens__hermann	Approve	approved as engineering manager	20 Dec 2010 12:46:56 GMT

* Printed versions are For Reference Only *



DICOM V3.0 Conformance Statement MammoWorkstation 4.6.0





© Image Diagnost International GmbH

Reproduction of this document without prior written approval is prohibited.
Image diagnost reserves the right to revise this manual

Image Diagnost International GmbH
Balanstraße 57
81541 München, Germany



Released



Document Versions

Version	Date	Author	Changes
1.00	04-Nov-03	Peter Heinlein	First Draft
1.01	10-Dec-03	Peter Heinlein	Minor corrections
1.10	07-Apr-04	Peter Heinlein	Added DICOM Worklist FIND & MPPS
1.35	10-Sept-04	Peter Heinlein	Minor corrections
1.5	24-Jan-05	Peter Heinlein	Added DR, CT, MR, Ultrasound, GSPS Documents
3.0	17-Jun-05	Peter Heinlein	Added Verification SCU and SCP
3.2	01-May-05	Peter Heinlein	Added SR Templates
4.5.3	01-Jul-10	Robert Lopuszanski	Added annex for jpeg compression Added encapsulated PDF Added storage commit push model
	Dec-2010	Siegfried Janssen	Updated Print
4.6.0	Dec-2010	Robert Lopuszanski	Reviewed and updated for 4.6.0
4.6.0	Dec-2010	Anca Dima	Measurement chapter reviewed and updated for 4.6.0



Table of Contents

0	INTRODUCTION	10
0.1	DEFINITION, ACRONYMS, ABBREVIATIONS	10
1	IMPLEMENTATION MODEL.....	12
1.1	APPLICATION DATA FLOW DIAGRAM	12
1.2	FUNCTIONAL DEFINITION OF APPLICATION ENTITIES (AE)	14
1.2.1	Verification AE	14
1.2.2	Worklist AE	14
1.2.3	Network AE.....	15
1.2.4	Print AE.....	16
1.2.5	Media Interchange AE	16
1.3	SEQUENCING OF REAL-WORLD ACTIVITIES	16
1.3.1	Scheduled Modality Performed Procedure Step	16
1.3.2	Manual Image Print.....	17
2	AE SPECIFICATIONS.....	19
2.1	VERIFICATION AE	19
2.1.1	Association establishment policies	19
2.1.1.1	General	19
2.1.1.2	Number of associations	19
2.1.1.3	Asynchronous nature.....	19
2.1.1.4	Implementation identifying information	19
2.1.2	Association initiation by real-world activity	19
2.1.2.1	Real-world activity: "Verification"	19
2.1.2.1.1	Associated real-world activity for "Request verification".....	19
2.1.2.1.2	Proposed presentation contexts for "Request verification"	20
2.1.3	Association Acceptance Policy.....	20
2.2	WORKLIST AE.....	20
2.2.1	Association establishment policies	20
2.2.1.1	General	20
2.2.1.2	Number of associations	20
2.2.1.3	Asynchronous nature.....	21
2.2.1.4	Implementation identifying information	21
2.2.2	Association initiation by real-world activity	21
2.2.2.1	Real-world activity: "Retrieve Scheduled Patient information"	21



- 2.2.2.1.1 Associated real-world activity for “Retrieve Scheduled Patient information” 21
- 2.2.2.1.2 Proposed presentation contexts for “Retrieve Scheduled Patient information” 21
- 2.2.2.2 Real-world activity: “Commit Worklist item” 23
- 2.2.2.2.1 Associated real-world activity for “Commit Worklist item” 23
- 2.2.2.2.2 Proposed presentation contexts for “Commit Worklist item” 23
- 2.2.3 Association Acceptance Policy 23
- 2.3 NETWORK AE 24
- 2.3.1 Association establishment policies 25
- 2.3.1.1 General 25
- 2.3.1.2 Number of associations 26
- 2.3.1.3 Asynchronous nature 26
- 2.3.1.4 Implementation identifying information 26
- 2.3.2 Association initiation by real-world activity 26
- 2.3.2.1 Real-world activity: “Query/Retrieve” 26
- 2.3.2.1.1 Associated real-world activity for “Query/Retrieve” 26
- 2.3.2.1.2 Association Initiation Policy 26
- 2.3.2.1.3 Proposed presentation contexts for “Query/Retrieve” 27
- 2.3.2.2 Real-world activity: “Send Images and Reports” 28
- 2.3.2.2.1 Associated real-world activity for “Send Images and Reports” 28
- 2.3.2.2.2 Proposed presentation contexts for “Send Images and Reports” 29
- 2.3.2.2.3 Storage commitment 32
- 2.3.3 Association Acceptance Policy 33
- 2.3.3.1 Real-world activity: “Receive Images and Documents” 33
- 2.3.3.1.1 Associated real-world activity for “Receive Images and Documents” 33
- 2.3.3.1.2 Accepted presentation contexts for “Receive Images and Documents” 33
- 2.3.3.1.3 VOI LUT Attributes Module 42
- 2.3.3.1.4 Presentation Context Acceptance Criterion 42
- 2.3.3.1.5 Storage commitment (storage commitment broker) 42
- 2.4 PRINT AE 43
- 2.4.1 Association establishment policies 43
- 2.4.1.1 General 43
- 2.4.1.2 Number of associations 43
- 2.4.1.3 Asynchronous nature 44
- 2.4.1.4 Implementation identifying information 44
- 2.4.2 Association initiation by real-world activity 44
- 2.4.2.1 Real-world activity: “Manual Image Print” 44
- 2.4.2.1.1 Associated real-world activity for “Manual Image Print” 44
- 2.4.2.1.2 Proposed Presentation contexts for “Manual Image Print” 44
- 2.4.3 Association Acceptance Policy 45



- 2.5 MEDIA INTERCHANGE AE46
 - 2.5.1 Introduction46
 - 2.5.2 Implementation Model46
 - 2.5.2.1 Application Data Flow Diagram46
 - 2.5.3 Functional definition of the Media Interchange AE.....46
 - 2.5.4 File Meta Information For Implementation Class And Version46
 - 2.5.5 AE SPECIFICATION46
 - 2.5.5.1 Media Interchange AE Specification46
 - 2.5.5.2 File Meta Information for the Application Entity.....47
 - 2.5.5.3 Real world activities for the Media Interchange AE.....47
 - 2.5.5.3.1 Display Directory47
 - 2.5.5.3.2 Real world activity “CD-R burning”47
 - 2.5.5.3.3 Real world activity “CD-R/CDROM reading”48
 - 2.5.6 Augmented And Private Application Profiles49
 - 2.5.7 Extensions, Specializations And Privatizations Of SOP Classes And Transfer Syntaxes49
- 3 COMMUNICATION PROFILES50**
 - 3.1 SUPPORTED COMMUNICATION STACKS50
 - 3.2 TCP/IP STACK.....50
 - 3.2.1 Physical Media Support.....50
 - 3.2.2 Internet Protocol Version 6 (IPv6)50
- 4 EXTENSIONS/SPECIALIZATIONS/PRIVATIZATIONS.....50**
 - 4.1 STANDARD EXTENDED/SPECIALIZED/PRIVATE SOPs.....50
 - 4.2 PRIVATE TRANSFER SYNTAXES50
- 5 CONFIGURATION50**
 - 5.1 AE TITLE/PRESENTATION ADDRESS MAPPING51
- 6 SUPPORT OF EXTENDED CHARACTER SETS51**
- 7 ANNEX: STRUCTURED REPORTING SOP CLASS DEFINITION52**
 - 7.1 INTRODUCTION52
 - 7.2 STRUCTURED REPORT SOP CLASS DEFINITIONS.....52
 - 7.2.1 Templates52
 - 7.2.1.1 TID 4200 BREAST IMAGING REPORT52
 - 7.2.1.2 TID 1204 LANGUAGE OF CONTENT ITEM AND DESCENDANTS.....52
 - 7.2.1.3 TID 4202 BREAST IMAGING REPORT NARRATIVE53
 - 7.2.1.4 TID 4208 BREAST IMAGING REPORT SUPPLEMENTARY DATA.....53
 - 7.2.1.5 TID 4205 BREAST COMPOSITION SECTION.....54



7.2.1.6 TID 4206 BREAST IMAGING REPORT FINDING SECTION54

7.2.1.7 TID 4201 BREAST IMAGING PROCEDURE REPORTED.....56

7.2.1.8 TID 4203 BREAST IMAGING ASSESSMENT.....57

7.2.1.9 TID 1400 LINEAR MEASUREMENT58

8 ANNEX: PRINT MANAGEMENT SOP CLASS DEFINITION59

8.1 INTRODUCTION 59

8.2 PRINT MANAGEMENT SOP CLASS DEFINITIONS 59

8.2.1 Basic Film Session SOP Class.....59

8.2.1.1 IOD Description59

8.2.1.1.1 IOD Modules 59

8.2.1.1.2 Basic Film Session Presentation Module..... 59

8.2.1.1.3 Basic Film Session Relationship Module..... 60

8.2.1.2 DIMSE Service Group60

8.2.1.2.1 N-CREATE 60

8.2.2 Basic Film Box SOP Class.....61

8.2.2.1 IOD Description61

8.2.2.1.1 IOD Modules 61

8.2.2.1.2 Basic Film Box Presentation Module..... 62

8.2.2.1.3 Basic Film Box Relationship Module 63

8.2.2.2 DIMSE Service Group64

8.2.2.2.1 N-CREATE 64

8.2.2.2.2 N-DELETE..... 66

8.2.2.2.3 N-ACTION..... 66

8.2.3 Image Box SOP Class67

8.2.3.1 Basic Grayscale Image Box SOP Class.....67

8.2.3.1.1 IOD Description 67

8.2.3.1.2 DIMSE Service Group 68

8.2.4 Printer SOP Class.....70

8.2.4.1 IOD Description70

8.2.4.1.1 IOD Modules 70

8.2.4.1.2 Printer Module..... 70

8.2.4.2 DIMSE Service Group71

8.2.4.2.1 N-EVENT-REPORT 71

8.2.4.2.2 N-GET..... 72

8.2.5 Presentation LUT SOP Class.....72

8.2.5.1 IOD Description73

8.2.5.1.1 IOD Modules 73

8.2.5.1.2 Presentation LUT Module 73

8.2.5.2 DIMSE Service Group73



8.2.5.2.1 N-CREATE 73

9 ANNEX: JPEG DECOMPRESSION.....75

9.1 TRANSFER SYNTAX UID 1.2.840.10008.1.2.4.50, JPEG BASELINE (PROCESS 1) 75

9.2 TRANSFER SYNTAX UID 1.2.840.10008.1.2.4.51, JPEG EXTENDED (PROCESS 2 & 4) 75

9.3 TRANSFER SYNTAX UID 1.2.840.10008.1.2.4.70, JPEG LOSSLESS, NON-HIERARCHICAL, FIRST-ORDER
PREDICTION (PROCESS 14 [SELECTION VALUE 1]) 76

9.4 TRANSFER SYNTAX UID 1.2.840.10008.1.2.4.91, JPEG 2000 IMAGE COMPRESSION..... 76

9.5 TRANSFER SYNTAX UID 1.2.840.10008.1.2.4.90, JPEG 2000 IMAGE COMPRESSION (LOSSLESS ONLY) 76



0 Introduction

This is a conformance statement for the “MammoWorkstation” mammography workstation. MammoWorkstation is a diagnostic softcopy reading workstation for digital or digitized mammograms and related modalities. It is able to exchange DICOM data via network and DICOMDIR. The report generation module includes the creation of DICOM text SR, DICOM GSPS and DICOM encapsulated PDF. The print module supports printing including presentation LUTs.

0.1 Definition, Acronyms, Abbreviations

This document uses the following abbreviations.

AE	Application Entity
CAD	Computer Aided Detection
CD-R	Compact Disk Recordable
CSE	Customer Service Engineer
CR	Computed Radiography
CT	Computed Tomography
DHCP	Dynamic Host Configuration Protocol
DICOM	Digital Imaging and Communications in Medicine
DIMSE	DICOM Message Service Element
DIT	Directory Information Tree (LDAP)
DN	Distinguished Name (LDAP)
DNS	Domain Name System
DX	Digital X-ray
FSC	File-Set Creator
FSU	File-Set Updater
FSR	File-Set Reader
GSDF	Grayscale Standard Display Function
GSPS	Grayscale Softcopy Presentation State
HIS	Hospital Information System
IHE	Integrating the Healthcare Enterprise
IOD	Information Object Definition
IPv4	Internet Protocol version 4
IPv6	Internet Protocol version 6
ISO	International Organization for Standards
IO	Intra-oral X-ray
JPEG	Joint Photographic Experts Group
LUT	Look-up Table
MAR	Medication Administration Record
MPEG	Moving Picture Experts Group
MG	Mammography (X-ray)
MPPS	Modality Performed Procedure Step
MR	Magnetic Resonance Imaging
MSPS	Modality Scheduled Procedure Step
MTU	Maximum Transmission Unit (IP)
MWL	Modality Worklist
NTP	Network Time Protocol
O	Optional (Key Attribute)
OSI	Open Systems Interconnection
PACS	Picture Archiving and Communication System
PDU	Protocol Data Unit
R	Required (Key Attribute)
RDN	Relative Distinguished Name (LDAP)



RIS	Radiology Information System
SC	Secondary Capture
SCP	Service Class Provider
SCU	Service Class User
SOP	Service-Object Pair
SPS	Scheduled Procedure Step
SR	Structured Reporting
TCP/IP	Transmission Control Protocol/Internet Protocol
U	Unique (Key Attribute)
UID	Unique Identifier
UL	Upper Layer
US	Ultrasound
VL	Visible Light
VR	Value Representation



1 Implementation model

1.1 Application data flow diagram

The application data flow models for the MammoWorkstation are shown in the following illustration:

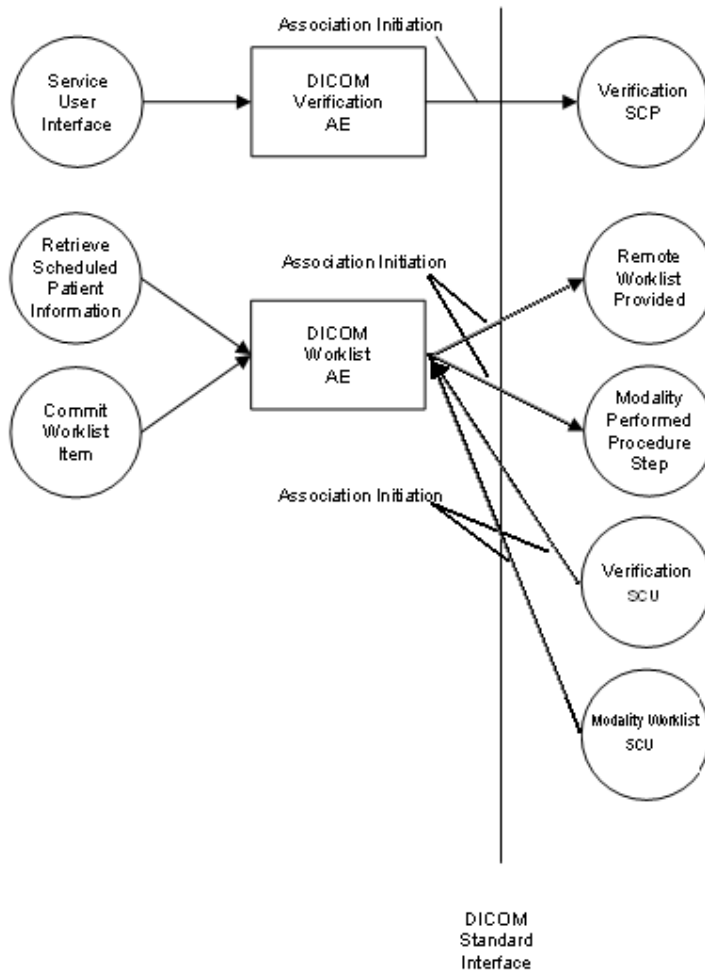


Figure 1.1-1 (Part1) – Application data flow diagram

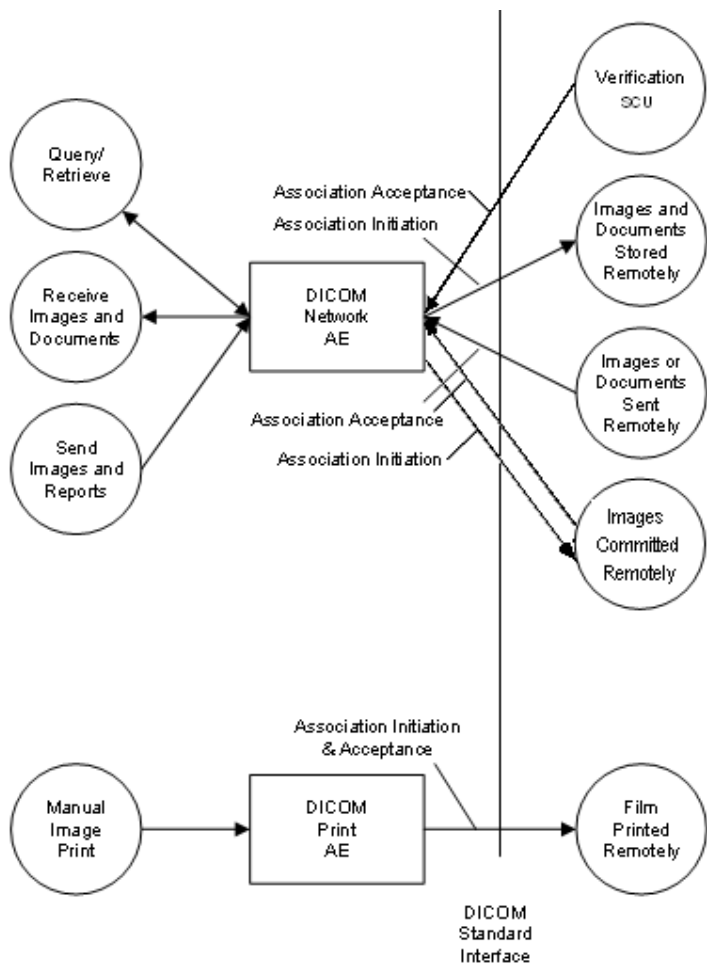


Figure 1.1-1 (Part2) – Application data flow diagram

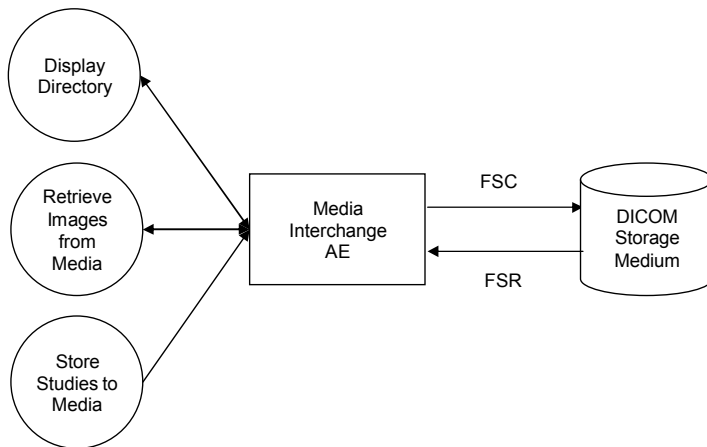




Figure 1.1-1 (Part3) – Application data flow diagram

The MammoWorkstation application is composed of the following Application Entities (AE):

1. DICOM Worklist Application Entity
2. DICOM Network AE
3. DICOM Print AE
4. DICOM Media Interchange AE

See chapter "Functional definition of Application Entities (AE)" for a detailed description.

1.2 Functional definition of Application Entities (AE)

1.2.1 Verification AE

The **DICOM Verification** Application Entity provides DICOM protocol verification of remote DICOM Services. The DICOM Verification AE is activated, when an user with administrative privileges opens the configuration settings user interface and activates the "Echo" function for a selected DICOM partner configuration.

The DICOM Verification AE provides the following Real World Activities:

- Verify the ability of a foreign DICOM application on a remote node to respond to DICOM messages (DICOM Verification SCU).

1.2.2 Worklist AE

The **DICOM Worklist** Application Entity provides DICOM protocol communication with remote DICOM Modality Worklist SCPs and SCUs.

The DICOM Worklist AE acting as SCU is activated when the user opens the worklist tab of the MammoWorkstation. He can search in the scheduled procedure steps inserted in a remote RIS using different filters (e.g. scheduled station, scheduled date). Then the scheduled procedure steps are shown in a list control of the MammoWorkstation GUI. The user can select one item and start to review the corresponding patient. The Performed Procedure Step Manager is used to update status information. The DICOM Worklist AE provides the following Real World Activities:

- Retrieve scheduled patient information from a remote Worklist Server (DICOM Modality Worklist SCU).
- Set Modality Performed Procedure Step

The DICOM Worklist AE acting as SCP for other SCUs is activated automatically when MammoWorkstation is started.

- Provide scheduled patient information from a remote worklist client (DICOM Modality Worklist SCP).
- Receive Modality Performed Procedure Step



1.2.3 Network AE

The **DICOM Network** Application Entity provides DICOM protocol communication for images, structured reports, GSPS data, encapsulated PDF and mammography CAD documents. The DICOM Network AE is automatically brought up when the MammoWorkstation is powered on.

The DICOM Network AE provides the following Real World Activities:

- Manual and automatic Query/Retrieve a remote DICOM archive to obtain a list of data at Patient level by selecting the remote DICOM AE (DICOM Query/Retrieve SCU). From the list the user can select an entry and retrieve the SOP Classes supported by the MammoWorkstation from the remote DICOM AE. The query is selective based on criteria described below in the document.
- Receive images, structured reports, GSPS data, encapsulated PDF and mammography CAD documents sent from a remote DICOM AE to the MammoWorkstation (DICOM Receive SCP).
- Send images, structured reports, GSPS data, encapsulated PDF and mammography CAD documents to remote DICOM Storage stations (DICOM Store SCU).
- Initiate storage commit for documents created on the workstation (Storage Commitment N-ACTION).
- Receive storage commit (Storage Commitment N-EVENT-REPORT)

The DICOM Network AE initiates the following operations:

- Access patient demographics and pixel data in the local database.
- Initiate a DICOM association to ask for remote patient demographics.
- Initiate storage commit for documents created on the workstation (Storage Commitment N-ACTION).
- Initiate a DICOM association to ask for transmit SOP Classes (images, structured reports, mammography CAD documents) from a remote DICOM AE to MammoWorkstation
- Initiate a DICOM association to send DICOM SOP Classes (secondary capture images, structured reports) to a remote DICOM AE.

The DICOM Network AE waits for association requests from Remote AE:

- Answer to DICOM associations transmitting DICOM SOP Classes (images, structured reports, mammography CAD documents) to be stored on the MammoWorkstation
- Receive storage commit (Storage Commitment N-EVENT-REPORT). For documents created on the workstation, the status is updated in the local image database. For documents sent to the MammoWorkstation, the storage commit answer is forwarded to the originator of the document.
- Receive storage commit request (Storage Commitment N-ACTION) for DICOM documents sent to the MammoWorkstation and forward the request to a configurable PACS archive. The MammoWorkstation is not an image archive, it does not answer storage commit requests directly.



1.2.4 Print AE

The **DICOM Print** Application Entity provides DICOM protocol communication with Remote DICOM Printers. The DICOM Print AE is activated when the user requests for a print.

The DICOM Print AE provides the following Real World Activities:

- Manual Image Print (DICOM Print SCU). For this operation the operator selects an image in the Viewer, then prints the image using the “Print” function.

The DICOM Print AE initiates the following operations:

- Initiate a DICOM association to send DICOM SOP Classes (corresponding to the DICOM Print Management service class) to a remote DICOM Printer.

1.2.5 Media Interchange AE

The **DICOM Media Interchange** Application Entity is activated when the user requests to browse a DICOM Media. The DICOM Media Interchange AE provides the following Real World Activities:

- Display Directory of a DICOM Storage Medium.
- Retrieve Images from a DICOM Storage Medium.
- Store studies to a DICOM Storage Medium.

Possible Media are local disk, CD or DVD.


The DICOM **Media Interchange AE** initiates the following operations:

- When the user clicks on a configured “DICOM CD” source, referring to an existing media, the user can browse the media and he can select and import studies from the media.
- When the user clicks on the “Export to DICOM CD” button, the user can create a new media containing the selected studies.

1.3 Sequencing of real-world activities

1.3.1 Scheduled Modality Performed Procedure Step

Sequencing of worklist, storage and performed procedure step management activities:

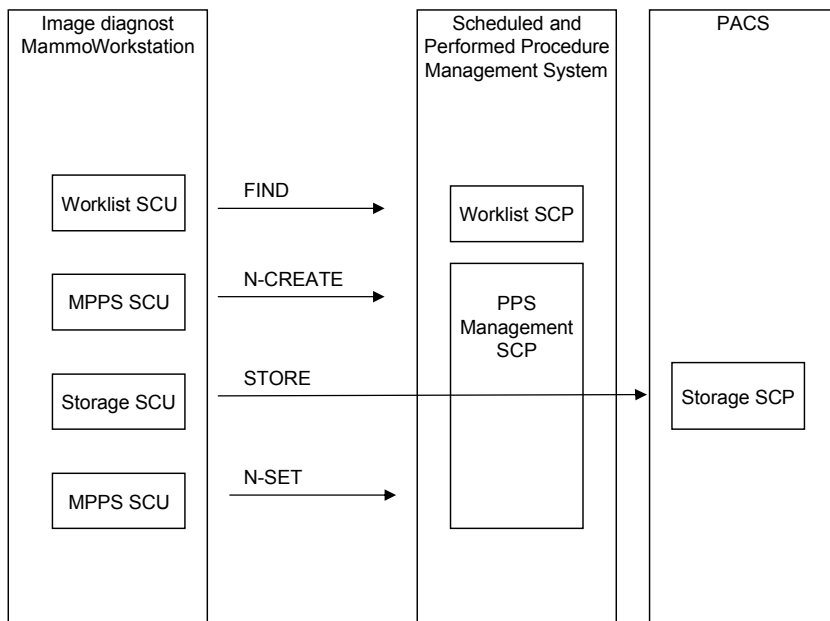

1.3.2 Manual Image Print

1. The user selects an image in the viewer and opens the Print Manager dialog.
2. The Print Manager dialog is displayed with the "Print" button initially locked. The Print manager retrieves the printer status by activating the DICOM Print SCU AE that initiates the following actions.
3. Initiates a DICOM association and selects a Presentation Context
4. N-GETs printer status from the Printer SOP Instance
 - a. If Printer Status is NORMAL
 - i. Association is released
 - ii. The Print manager enables the "Print" button, i.e. the user may initiate printing.
 - b. Else
 - i. Signal printer failure to the user
 - ii. The "Print" function remains locked
 - iii. Association is aborted
5. The user may select options in the Print Manager such as applying a calibration LUT and then prints the image using the "Print" button.
6. The Print Manager receives the "Print" request, composes a film, and activates the DICOM Print AE again. The AE then initiates the following actions:
7. Initiate a DICOM association and select the Presentation Context
8. N-CREATEs a Basic Film Session SOP instance
9. N-CREATEs a Basic Film Box SOP instance for the current film
10. N-SETs the Basic Film Box SOP Instance with the Image Box SOP Instance for each image on the film.
11. If the user selected a calibration LUT:
 - a. composes the corresponding Presentation LUT
 - b. N-CREATEs a Presentation LUT SOP instance for the current film
 - c. Creates a Referenced Presentation LUT Sequence in the Basic Film SOP Instance that contains a reference to this Presentation LUT SOP instance



12. N-ACTIONS on the Basic Film Box SOP Instance
13. N-DELETES on the Basic Film Box SOP Instance
14. Releases the DICOM association after printing is successful or failure has been signaled to the user.



2 AE specifications

2.1 Verification AE

The Verification AE, provides Standard Conformance to the following DICOM V3.0 Service Object Pair (SOP) Class as a Verification Service Class User (SCU):

SOP Class Name	SOP Class UID
Verification	1.2.840.10008.1.1

2.1.1 Association establishment policies

2.1.1.1 General

The VERIFICATION AE will initiate an association for DIMSE C-ECHO service operations.

2.1.1.2 Number of associations

The VERIFICATION AE allows a single association for association initiation.

2.1.1.3 Asynchronous nature

The VERIFICATION AE does not support asynchronous communication (multiple outstanding transactions over a single association).

2.1.1.4 Implementation identifying information

The Implementation Class Unique Identifier (UID) for the VERIFICATION AE is:

Implementation UID	1.2.276.0.69.25.1.1
--------------------	---------------------

2.1.2 Association initiation by real-world activity

The following real world activities initiate association for the VERIFICATION AE:

1. Verify the ability of a foreign DICOM application on a remote node to respond to DICOM messages

2.1.2.1 Real-world activity: "Verification"

2.1.2.1.1 Associated real-world activity for "Request verification"

The DICOM Verification AE is activated when an user with administrative privileges opens the configuration settings user interface and activates the "Echo" function for selected DICOM partner configuration. The status of the C-ECHO Response from the remote application is indicated in the user interface. A GUI control indicates "success" and all other status.



2.1.2.1.2 Proposed presentation contexts for “Request verification”

The presentation contexts that are proposed by VERIFICATION AE for the verification operation are:

Presentation Context Table					
Abstract Syntax		Transfer Syntax		Role	Extended Negotiation
Name	UID	Name List	UID List		
Verification	1.2.840.10008.1.1	Explicit VR Little Endian	1.2.840.10008.1.2.1	SCU	None
		Implicit VR Little Endian	1.2.840.10008.1.2		
		Explicit VR Big Endian	1.2.840.10008.1.2.2		

2.1.3 Association Acceptance Policy

VERIFICATION AE does not respond to attempts by a remote AE to open an association.

2.2 Worklist AE

The Worklist AE provides Standard Conformance to the following DICOM V3.0 Service Object Pair (SOP) Class as a Modality Worklist Service Class User (SCU):

SOP Class Name	SOP Class UID
Verification	1.2.840.10008.1.1
Modality Worklist Information Model – FIND	1.2.840.10008.5.1.4.31
Modality Performed Procedure Step SOP Class	1.2.840.10008.3.1.2.3.3

2.2.1 Association establishment policies

2.2.1.1 General

The WORKLIST AE will initiate an association as a Modality Worklist Service Class User requesting Application Entity Title and patient data.

The PDU size is 64234 bytes.

2.2.1.2 Number of associations

The WORKLIST AE opens one association for each request. The WORKLIST AE will accept up to ten associations at the same time



2.2.1.3 Asynchronous nature

The WORKLIST AE does not support asynchronous communication (multiple outstanding transactions over a single association).

2.2.1.4 Implementation identifying information

The Implementation Class Unique Identifier (UID) for the WORKLIST AE is:

Implementation UID	1.2.276.0.69.25.1.1
--------------------	---------------------

2.2.2 Association initiation by real-world activity

The following real world activities initiate association for the WORKLIST AE:

1. Retrieve Scheduled Patient information
2. User Commits Worklist item

2.2.2.1 Real-world activity: “Retrieve Scheduled Patient information”

2.2.2.1.1 Associated real-world activity for “Retrieve Scheduled Patient information”

When the user clicks on the “Worklist” tab from the MammoWorkstation and then on the “Refresh” button, the Worklist AE initiates an association for the appropriate Modality Worklist Service Class that corresponds to the set of data requested to be transferred. Then the Worklist AE performs a C-FIND and it is closed when all the response items are received from the remote SCP.

2.2.2.1.2 Proposed presentation contexts for “Retrieve Scheduled Patient information”

The presentation contexts that are proposed by WORKLIST AE for the FIND operation are:

Presentation Context Table					
Abstract Syntax		Transfer Syntax		Role	Extended Negotiation
Name	UID	Name List	UID List		
Modality Worklist Find	1.2.840.10008.5.1.4.31	Explicit VR Little Endian	1.2.840.10008.1.2.1	SCU	None
		Implicit VR Little Endian	1.2.840.10008.1.2		
		Explicit VR Big Endian	1.2.840.10008.1.2.2		

2.2.2.1.2.1 SOP Specific Conformance Statement for Modality Worklist Find



WORKLIST AE provides standard conformance as an SCU to the Modality Worklist Service Class and uses the following elements for this SOP class.:

Description / Name	Tag	Matching Key Type	Return Key Type
SOP Common			
Specific Character Set	(0008,0005)	O	1C
Scheduled Procedure Step			
Scheduled Procedure Step Sequence	(0040,0100)	R	1
>Scheduled Station AE Title	(0040,0001)	R	1
>Scheduled Procedure Step Start Date	(0040,0002)	R	1
>Scheduled Procedure Step Start Time	(0040,0002)	R	1
>Modality	(0008,0060)	R	1
>Scheduled Performing Physician's Name	(0040,0006)	R	2
>Scheduled Procedure Step Description	(0040,0007)	O	1C
>Scheduled Station Name	(0040,0010)	O	2
>Scheduled Procedure Step ID	(0040,0009)	O	1
Requested Procedure			
Requested Procedure ID	(0040,1001)	O	1
Requested Procedure Description	(0032,1060)	O	1C
Study Instance UID	(0020,000D)	O	1
Imaging Service Request			
Accession Number	(0008,0050)	O	2
Referring Physician's Name	(0008,0090)	O	2
Patient Identification			



Patient's Name	(0010,0010)	R	1
Patient ID	(0010,0020)	R	1
Patient Demographic			
Patient's Birth Date	(0010,0030)	O	2
Patient's Sex	(0010,0040)	O	2

2.2.2.2 Real-world activity: "Commit Worklist item"

2.2.2.2.1 Associated real-world activity for "Commit Worklist item"

When the user select an item from the worklist, the MammoWorkstation sends a N-CREATE request.

2.2.2.2.2 Proposed presentation contexts for "Commit Worklist item"

The presentation contexts that are proposed by WORKLIST AE for the FIND operation are:

Presentation Context Table					
Abstract Syntax		Transfer Syntax		Role	Extended Negotiation
Name	UID	Name List	UID List		
Modality Performed Procedure Step	1.2.840.10008.3.1.2.3.3	Explicit VR Little Endian	1.2.840.10008.1.2.1	SCU	None
		Implicit VR Little Endian	1.2.840.10008.1.2		
		Explicit VR Big Endian	1.2.840.10008.1.2.2		

2.2.3 Association Acceptance Policy

When the WORKLIST AE accepts an association, it will process modality worklist find and modality performed procedure step request.



2.3 NETWORK AE

The Network AE provides Standard Conformance to the following DICOM V3.0 Service Object Pair (SOP) Classes as a Service Class User (SCU):

SOP Class Name	SOP Class UID
Verification	1.2.840.10008.1.1
Computed Radiography Image Storage	1.2.840.10008.5.1.4.1.1.1
Digital X-Ray Image Storage – For Presentation	1.2.840.10008.5.1.4.1.1.1.1
Digital X-Ray Image Storage – For Processing	1.2.840.10008.5.1.4.1.1.1.1.1
Digital Mammography X-Ray Image Storage – For Presentation	1.2.840.10008.5.1.4.1.1.1.2
Digital Mammography X-Ray Image Storage – For Processing	1.2.840.10008.5.1.4.1.1.1.2.1
CT Image Storage	1.2.840.10008.5.1.4.1.1.2
MR Image Storage	1.2.840.10008.5.1.4.1.1.4
US Image Storage	1.2.840.10008.5.1.4.1.1.6.1
Secondary Capture Image Storage	1.2.840.10008.5.1.4.1.1.7
Grayscale Softcopy Presentation State Storage	1.2.840.10008.5.1.4.1.1.11.1
Basic Text SR	1.2.840.10008.5.1.4.1.1.88.11
Enhanced Structured Reporting	1.2.840.10008.5.1.4.1.1.88.22
Comprehensive Structured Reporting	1.2.840.10008.5.1.4.1.1.88.33
Mammography CAD SR	1.2.840.10008.5.1.4.1.1.88.50
Patient Root Query/Retrieve Information Model – FIND	1.2.840.10008.5.1.4.1.2.1.1
Patient Root Query/Retrieve Information Model - MOVE	1.2.840.10008.5.1.4.1.2.1.2
Study Root Query/Retrieve Information Model – FIND	1.2.840.10008.5.1.4.1.2.2.1
Study Root Query/Retrieve Information Model - MOVE	1.2.840.10008.5.1.4.1.2.2.2
Encapsulated PDF Storage	1.2.840.10008.5.1.4.1.1.104.1



Storage Commitment Push Model	1.2.840.10008.1.20.1
-------------------------------	----------------------

The Network AE provides Standard Conformance to the following DICOM V3.0 Service Object Pair (SOP) Classes as a Service Class Provider (SCP):

SOP Class Name	SOP Class UID
Verification	1.2.840.10008.1.1
Computed Radiography Image Storage	1.2.840.10008.5.1.4.1.1.1
Digital X-Ray Image Storage – For Presentation	1.2.840.10008.5.1.4.1.1.1.1
Digital X-Ray Image Storage – For Processing	1.2.840.10008.5.1.4.1.1.1.1.1
Digital Mammography X-Ray Image Storage – For Presentation	1.2.840.10008.5.1.4.1.1.1.2
Digital Mammography X-Ray Image Storage – For Processing	1.2.840.10008.5.1.4.1.1.1.2.1
CT Image Storage	1.2.840.10008.5.1.4.1.1.2
MR Image Storage	1.2.840.10008.5.1.4.1.1.4
US Image Storage	1.2.840.10008.5.1.4.1.1.6.1
Secondary Capture Image Storage	1.2.840.10008.5.1.4.1.1.7
Grayscale Softcopy Presentation State Storage	1.2.840.10008.5.1.4.1.1.11.1
Basic Text SR	1.2.840.10008.5.1.4.1.1.88.11
Enhanced Structured Reporting	1.2.840.10008.5.1.4.1.1.88.22
Comprehensive Structured Reporting	1.2.840.10008.5.1.4.1.1.88.33
Mammography CAD SR	1.2.840.10008.5.1.4.1.1.88.50
Encapsulated PDF Storage	1.2.840.10008.5.1.4.1.1.104.1
Storage Commitment Push Model	1.2.840.10008.1.20.1

2.3.1 Association establishment policies

2.3.1.1 General



The Network AE will initiate an association as an SCU of Storage Services when

- a local operator requests to send images or structured reports over the network to a remote Storage Service Class provider.

The PDU size is 64234bytes.

2.3.1.2 Number of associations

The Network AE opens one Store association for each image to be sent to the storage station. The Network AE will accept up to five associations at the same time.

2.3.1.3 Asynchronous nature

The Network AE does not support asynchronous communication (multiple outstanding transactions over a single association).

2.3.1.4 Implementation identifying information

The Implementation Class Unique Identifier (UID) for the NETWORK AE is:

Implementation UID	1.2.276.0.69.25.1.1
--------------------	---------------------

2.3.2 Association initiation by real-world activity

2.3.2.1 Real-world activity: “Query/Retrieve”

2.3.2.1.1 Associated real-world activity for “Query/Retrieve”

The operator queries a Remote database by selecting the “search” tab and clicking on the icon representing the DICOM Remote AE. A list control displays the results upon successful query.

Then the operator can select one or several Patients and can either click the “Prefetch” icon to retrieve the selection on the MammoWorkstation database or click the “Select” icon to retrieve the selection on the MammoWorkstation and start reviewing the first Patient in the list.

These operations will cause:

- The DICOM NETWORK AE to initiate a DICOM association
- The DICOM NETWORK AE to emit a C-FIND request to get a list of patients regarding the criteria listed below, to get the appropriate studies, series or images for the selected list of patients.
- The DICOM NETWORK AE to emit a C-FIND request to specify a selected list of studies/series/images to be sent by the Remote Host to the MammoWorkstation

2.3.2.1.2 Association Initiation Policy

The NETWORK AE initiates a new association for querying Patient Folders (or Studies/Series/Images) on a remote DICOM AE. This association corresponds to one Real World Activity:



2.3.2.1.3 Proposed presentation contexts for “Query/Retrieve”

The presentation contexts that are proposed by NETWORK AE for the FIND and MOVE operation are:

Presentation Context Table					
Abstract Syntax		Transfer Syntax		Role	Extended Negotiation
Name	UID	Name List	UID List		
Patient Root Query/Retrieve Information Model - FIND	1.2.840.10008.5.1.4.1.1.1	Explicit VR Little Endian	1.2.840.10008.1.2.1	SCU	None
		Implicit VR Little Endian	1.2.840.10008.1.2		
		Explicit VR Big Endian	1.2.840.10008.1.2.2		
Patient Root Query/Retrieve Information Model - MOVE	1.2.840.10008.5.1.4.1.1.2	Explicit VR Little Endian	1.2.840.10008.1.2.1	SCU	None
		Implicit VR Little Endian	1.2.840.10008.1.2		
		Explicit VR Big Endian	1.2.840.10008.1.2.2		
Study Root Query/Retrieve Information Model – FIND	1.2.840.10008.5.1.4.1.2.1	Explicit VR Little Endian	1.2.840.10008.1.2.1	SCU	None
		Implicit VR Little Endian	1.2.840.10008.1.2		
		Explicit VR Big Endian	1.2.840.10008.1.2.2		
Study Root Query/Retrieve Information Model – MOVE	1.2.840.10008.5.1.4.1.2.2	Explicit VR Little Endian	1.2.840.10008.1.2.1	SCU	None
		Implicit VR Little Endian	1.2.840.10008.1.2		
		Explicit VR Big Endian	1.2.840.10008.1.2.2		

2.3.2.1.3.1 SOP Specific Conformance Statement for the Model Patient Root Query/Retrieve Information Model – FIND SOP Class

Following are the status codes that are more specifically processed when receiving messages from a Query SCP equipment:

Service Status	Status Codes	Further Meaning	Application Behaviour When receiving Status Codes	Related Fields Processed if received
Refused	A700	Out of resources	Association is closed. Appropriate message is displayed to the user.	(0000,0902)



	0122	SOP Class not Supported	Association is closed. Appropriate message is displayed to the user.	(0000,0902)
Failed	A900	Identifier does not match SOP Class	Association is closed. Error message is displayed to the user.	(0000,0901) (0000,0902)
	Cxxx	Unable to process	Association is closed. Error message is displayed to the user.	(0000,0901) (0000,0902)
Cancel	FE00	Matching terminated due to cancel	Association is closed. Error message is displayed to the user.	None
Success	0000	Matching is complete – No final identifier is supplied		None
Pending	FF00	Matches are continuing – Current Match is supplied and any Optional Keys were supported in the same manner as Required Keys		Identifier
	FF01	Matches are continuing – Warning that one or more Optional Keys were not supported for existence and/or matching for this Identifier.		Identifier

2.3.2.1.3.2 SOP Specific Conformance Statement for the Model Study Root Query/Retrieve Information Model – FIND SOP Class

See section on Patient Root Query/Retrieve Information Model – FIND SOP Class.

2.3.2.2 Real-world activity: “Send Images and Reports”

2.3.2.2.1 Associated real-world activity for “Send Images and Reports”

When

1. the user selects some studies/series/images, clicks on the “Send” button from the MammoWorkstation database and then selects the remote AE,
2. the user closes a case with report data or annotations entered



- 3. the user creates a snapshot

the NETWORK AE initiates an association for each image, report, annotation selected or created. Every association is closed when the DICOM data has been sent (successfully or not) to the remote SCP. The NETWORK AE supports lossless and lossy JPEG 2000 compression. Compression can be configured in the MammoWorkstation and is invoked automatically for images.

2.3.2.2.2 Proposed presentation contexts for “Send Images and Reports”

The presentation contexts that are proposed by NETWORK AE for the “Manual Storing” real world activity are:

Presentation Context Table - Proposed					
Abstract Syntax		Transfer Syntax		Role	Extended Negotiation
Name	UID	Name List	UID List		
Verification	1.2.840.10008.1	Implicit VR Little Endian	1.2.840.10008.1.1	SCU	None
Computed Radiography Image Storage	1.2.840.10008.5.1.4.1.1.1	Explicit VR Little Endian	1.2.840.10008.1.2.1	SCU	None
		Implicit VR Little Endian	1.2.840.10008.1.2		
		Explicit VR Big Endian	1.2.840.10008.1.2.2		
		JPEG 2000 Image Compression (Lossless Only)	1.2.840.10008.1.2.4.90		
		JPEG 2000 Image Compression	1.2.840.10008.1.2.4.91		
Digital X-Ray Image Storage – For Presentation	1.2.840.10008.5.1.4.1.1.1.1	Explicit VR Little Endian	1.2.840.10008.1.2.1	SCU	None
		Implicit VR Little Endian	1.2.840.10008.1.2		
		Explicit VR Big Endian	1.2.840.10008.1.2.2		
		JPEG 2000 Image Compression (Lossless Only)	1.2.840.10008.1.2.4.90		
		JPEG 2000 Image Compression	1.2.840.10008.1.2.4.91		
Digital X-Ray Image Storage – For Processing	1.2.840.10008.5.1.4.1.1.1.1.1	Explicit VR Little Endian	1.2.840.10008.1.2.1	SCU	None
		Implicit VR Little Endian	1.2.840.10008.1.2		
		Explicit VR Big Endian	1.2.840.10008.1.2.2		



		JPEG 2000 Image Compression (Lossless Only)	1.2.840.10008.1.2.4.90		
		JPEG 2000 Image Compression	1.2.840.10008.1.2.4.91		
Digital Mammography X-Ray Image Storage – For Presentation	1.2.840.10008.5.1.4.1.1.1.2	Explicit VR Little Endian	1.2.840.10008.1.2.1	SCU	None
		Implicit VR Little Endian	1.2.840.10008.1.2		
		Explicit VR Big Endian	1.2.840.10008.1.2.2		
		JPEG 2000 Image Compression (Lossless Only)	1.2.840.10008.1.2.4.90		
		JPEG 2000 Image Compression	1.2.840.10008.1.2.4.91		
Digital Mammography X-Ray Image Storage – For Processing	1.2.840.10008.5.1.4.1.1.1.2.1	Explicit VR Little Endian	1.2.840.10008.1.2.1	SCU	None
		Implicit VR Little Endian	1.2.840.10008.1.2		
		Explicit VR Big Endian	1.2.840.10008.1.2.2		
		JPEG 2000 Image Compression (Lossless Only)	1.2.840.10008.1.2.4.90		
		JPEG 2000 Image Compression	1.2.840.10008.1.2.4.91		
CT Image Storage	1.2.840.10008.5.1.4.1.1.2	Explicit VR Little Endian	1.2.840.10008.1.2.1	SCU	None
		Implicit VR Little Endian	1.2.840.10008.1.2		
		Explicit VR Big Endian	1.2.840.10008.1.2.2		
		JPEG 2000 Image Compression (Lossless Only)	1.2.840.10008.1.2.4.90		
		JPEG 2000 Image Compression	1.2.840.10008.1.2.4.91		
MR Image Storage	1.2.840.10008.5.1.4.1.1.4	Explicit VR Little Endian	1.2.840.10008.1.2.1	SCU	None
		Implicit VR Little Endian	1.2.840.10008.1.2		



		Explicit VR Big Endian	1.2.840.10008.1.2.2		
		JPEG 2000 Image Compression (Lossless Only)	1.2.840.10008.1.2.4.90		
		JPEG 2000 Image Compression	1.2.840.10008.1.2.4.91		
US Image Storage	1.2.840.10008.5.1.4.1.1.6.1	Explicit VR Little Endian	1.2.840.10008.1.2.1	SCU	None
		Implicit VR Little Endian	1.2.840.10008.1.2		
		Explicit VR Big Endian	1.2.840.10008.1.2.2		
		JPEG 2000 Image Compression (Lossless Only)	1.2.840.10008.1.2.4.90		
		JPEG 2000 Image Compression	1.2.840.10008.1.2.4.91		
Secondary Capture Image Storage	1.2.840.10008.5.1.4.1.1.7	Explicit VR Little Endian	1.2.840.10008.1.2.1	SCU	None
		Implicit VR Little Endian	1.2.840.10008.1.2		
		Explicit VR Big Endian	1.2.840.10008.1.2.2		
		JPEG 2000 Image Compression (Lossless Only)	1.2.840.10008.1.2.4.90		
		JPEG 2000 Image Compression	1.2.840.10008.1.2.4.91		
Grayscale Softcopy Presentation State Storage	1.2.840.10008.5.1.4.1.1.1.1.1	Explicit VR Little Endian	1.2.840.10008.1.2.1	SCU	None
		Implicit VR Little Endian	1.2.840.10008.1.2		
		Explicit VR Big Endian	1.2.840.10008.1.2.2		
		JPEG 2000 Image Compression (Lossless Only)	1.2.840.10008.1.2.4.90		
		JPEG 2000 Image Compression	1.2.840.10008.1.2.4.91		
Basic Text SR	1.2.840.10008.5.1.4.1.1.8	Explicit VR Little Endian	1.2.840.10008.1.2.1	SCU	None



	8.11	Implicit VR Little Endian	1.2.840.10008.1.2		
		Explicit VR Big Endian	1.2.840.10008.1.2.2		
Mammography CAD SR	1.2.840.10008.5.1.4.1.1.8 8.50	Explicit VR Little Endian	1.2.840.10008.1.2.1	SCU	None
		Implicit VR Little Endian	1.2.840.10008.1.2		
		Explicit VR Big Endian	1.2.840.10008.1.2.2		
Encapsulated PDF Storage	1.2.840.10008.5.1.4. 1.1.104.1	Explicit VR Little Endian	1.2.840.10008.1.2.1	SCU	None
		Implicit VR Little Endian	1.2.840.10008.1.2		
		Explicit VR Big Endian	1.2.840.10008.1.2.2		

2.3.2.2.2.1 SOP Specific Conformance Statement for Digital Mammography Image Storage – For Presentation SOP Class

The NETWORK AE will react in the following ways to the different C-STORE response status:

- Refused: the SCU will show the FAILED status for the image to be sent.
- Error: the SCU will show the FAILED status for the image to be sent.
- Warning: the SCU will show the SUCCESS status for the image to be sent.
- Success: the SCU will show the SUCCESS status for the image to be sent.

NETWORK AE will not attempt any extended negotiation.

2.3.2.2.2.2 SOP Specific Conformance Statement for Digital Mammography Image Storage – For Processing SOP Class

Same as SOP Specific Conformance Statement for Digital Mammography Image Storage – For Presentation SOP Class.

2.3.2.2.3 Storage commitment

The Storage commitment is initiated only for documents created on the MammoWorkstation when the user closes a case with report data or annotations entered or when the user creates a snapshot. It is not initiated when the user selects some studies/series/images, clicks on the “Send” button from the MammoWorkstation database and then selects the remote AE.

The declaration of remote host supporting Storage Commitment is done through a specific menu (DICOM Partner menu). When a successful commitment for some documents is received, the NETWORK AE flags them as committed in the local database.

The NETWORK AE initiates the following operations:

- Initiate a DICOM association to ask a remote host (Storage Commitment SCP) on specific images.



The NETWORK AE waits for association requests from remote storage commitment AE:

- Answer to DICOM associations transmitting Storage Commitment Notification (N-EVENT-REPORT)

2.3.3 Association Acceptance Policy

When the NETWORK AE accepts an association, it will receive any DICOM data transmitted on that association and store the supported SOP Classes on disk. By default, any Remote DICOM AE can send images to the NETWORK AE.

2.3.3.1 Real-world activity: “Receive Images and Documents”

The NETWORK AE accepts an association when it receives a valid association request from a DICOM Storage SCU.

2.3.3.1.1 Associated real-world activity for “Receive Images and Documents”

The NETWORK AE waits for any association. No operator is required to receive an image. If images are received for the currently opened patient context, these images are opened automatically.

2.3.3.1.2 Accepted presentation contexts for “Receive Images and Documents”

Presentation Context Table – Accepted					
Abstract Syntax		Transfer Syntax		Role	Extended Negotiation
Name	UID	Name List	UID List		
Verification	1.2.840.10008.1	Implicit VR Little Endian	1.2.840.10008.1.1	SCP	None
Computed Radiography Image Storage	1.2.840.10008.5.1.4.1.1.1	Explicit VR Little Endian	1.2.840.10008.1.2.1	SCP	None
		Implicit VR Little Endian	1.2.840.10008.1.2		
		Explicit VR Big Endian	1.2.840.10008.1.2.2		
		JPEG 2000 Image Compression (Lossless Only)	1.2.840.10008.1.2.4.90		
		JPEG 2000 Image Compression	1.2.840.10008.1.2.4.91		
		JPEG Baseline Lossy 8-Bit	1.2.840.10008.1.2.4.50		
		JPEG Baseline Lossy 12-Bit	1.2.840.10008.1.2.4.51		
		JPEG Lossless Hierarch., First-order prediction	1.2.840.10008.1.2.4.70		
Digital X-Ray Image Storage – For Presentation	1.2.840.10008.5.1.4.1.1.1.1	Explicit VR Little Endian	1.2.840.10008.1.2.1	SCP	None
		Implicit VR Little Endian	1.2.840.10008.1.2		



		Explicit VR Big Endian	1.2.840.10008.1.2.2		
		JPEG 2000 Image Compression (Lossless Only)	1.2.840.10008.1.2.4.90		
		JPEG 2000 Image Compression	1.2.840.10008.1.2.4.91		
		JPEG Baseline Lossy 8-Bit	1.2.840.10008.1.2.4.50		
		JPEG Baseline Lossy 12-Bit	1.2.840.10008.1.2.4.51		
		JPEG Lossless Hierarch., First-order prediction	1.2.840.10008.1.2.4.70		
Digital X-Ray Image Storage – For Processing	1.2.840.10008.5.1.4.1.1.1.1.1	Explicit VR Little Endian	1.2.840.10008.1.2.1	SCP	None
		Implicit VR Little Endian	1.2.840.10008.1.2		
		Explicit VR Big Endian	1.2.840.10008.1.2.2		
		JPEG 2000 Image Compression (Lossless Only)	1.2.840.10008.1.2.4.90		
		JPEG 2000 Image Compression	1.2.840.10008.1.2.4.91		
		JPEG Baseline Lossy 8-Bit	1.2.840.10008.1.2.4.50		
		JPEG Baseline Lossy 12-Bit	1.2.840.10008.1.2.4.51		
		JPEG Lossless Hierarch., First-order prediction	1.2.840.10008.1.2.4.70		
Digital Mammography X-Ray Image Storage – For Presentation	1.2.840.10008.5.1.4.1.1.1.2	Explicit VR Little Endian	1.2.840.10008.1.2.1	SCP	None
		Implicit VR Little Endian	1.2.840.10008.1.2		
		Explicit VR Big Endian	1.2.840.10008.1.2.2		
		JPEG 2000 Image Compression (Lossless Only)	1.2.840.10008.1.2.4.90		
		JPEG 2000 Image Compression	1.2.840.10008.1.2.4.91		
		JPEG Baseline Lossy 8-Bit	1.2.840.10008.1.2.4.50		
		JPEG Baseline Lossy 12-Bit	1.2.840.10008.1.2.4.51		



		JPEG Lossless Hierarch., First-order prediction	1.2.840.10008.1.2.4.70		
Digital Mammography X-Ray Image Storage – For Processing	1.2.840.10008.5.1.4.1.1.1.2.1	Explicit VR Little Endian	1.2.840.10008.1.2.1	SCP	None
		Implicit VR Little Endian	1.2.840.10008.1.2		
		Explicit VR Big Endian	1.2.840.10008.1.2.2		
		JPEG 2000 Image Compression (Lossless Only)	1.2.840.10008.1.2.4.90		
		JPEG 2000 Image Compression	1.2.840.10008.1.2.4.91		
		JPEG Baseline Lossy 8-Bit	1.2.840.10008.1.2.4.50		
		JPEG Baseline Lossy 12-Bit	1.2.840.10008.1.2.4.51		
		JPEG Lossless Hierarch., First-order prediction	1.2.840.10008.1.2.4.70		
CT Image Storage	1.2.840.10008.5.1.4.1.1.2	Explicit VR Little Endian	1.2.840.10008.1.2.1	SCP	None
		Implicit VR Little Endian	1.2.840.10008.1.2		
		Explicit VR Big Endian	1.2.840.10008.1.2.2		
		JPEG 2000 Image Compression (Lossless Only)	1.2.840.10008.1.2.4.90		
		JPEG 2000 Image Compression	1.2.840.10008.1.2.4.91		
		JPEG Baseline Lossy 8-Bit	1.2.840.10008.1.2.4.50		
		JPEG Baseline Lossy 12-Bit	1.2.840.10008.1.2.4.51		
		JPEG Lossless Hierarch., First-order prediction	1.2.840.10008.1.2.4.70		
MR Image Storage	1.2.840.10008.5.1.4.1.1.4	Explicit VR Little Endian	1.2.840.10008.1.2.1	SCP	None
		Implicit VR Little Endian	1.2.840.10008.1.2		
		Explicit VR Big Endian	1.2.840.10008.1.2.2		
		JPEG 2000 Image Compression (Lossless Only)	1.2.840.10008.1.2.4.90		



		JPEG 2000 Image Compression	1.2.840.10008.1.2.4.91		
		JPEG Baseline Lossy 8-Bit	1.2.840.10008.1.2.4.50		
		JPEG Baseline Lossy 12-Bit	1.2.840.10008.1.2.4.51		
		JPEG Lossless Hierarch., First-order prediction	1.2.840.10008.1.2.4.70		
US Image Storage	1.2.840.10008.5.1.4.1.1.6.1	Explicit VR Little Endian	1.2.840.10008.1.2.1	SCP	None
		Implicit VR Little Endian	1.2.840.10008.1.2		
		Explicit VR Big Endian	1.2.840.10008.1.2.2		
		JPEG 2000 Image Compression (Lossless Only)	1.2.840.10008.1.2.4.90		
		JPEG 2000 Image Compression	1.2.840.10008.1.2.4.91		
		JPEG Baseline Lossy 8-Bit	1.2.840.10008.1.2.4.50		
		JPEG Baseline Lossy 12-Bit	1.2.840.10008.1.2.4.51		
		JPEG Lossless Hierarch., First-order prediction	1.2.840.10008.1.2.4.70		
Secondary Capture Image Storage	1.2.840.10008.5.1.4.1.1.7	Explicit VR Little Endian	1.2.840.10008.1.2.1	SCP	None
		Implicit VR Little Endian	1.2.840.10008.1.2		
		Explicit VR Big Endian	1.2.840.10008.1.2.2		
		JPEG 2000 Image Compression (Lossless Only)	1.2.840.10008.1.2.4.90		
		JPEG 2000 Image Compression	1.2.840.10008.1.2.4.91		
		JPEG Baseline Lossy 8-Bit	1.2.840.10008.1.2.4.50		
		JPEG Baseline Lossy 12-Bit	1.2.840.10008.1.2.4.51		
		JPEG Lossless Hierarch., First-order prediction	1.2.840.10008.1.2.4.70		
Grayscale Softcopy Presentation State Storage	1.2.840.10008.5.1.4.1.1.1.1.1	Explicit VR Little Endian	1.2.840.10008.1.2.1	SCP	None
		Implicit VR Little Endian	1.2.840.10008.1.2		



		Explicit VR Big Endian	1.2.840.10008.1.2.2		
		JPEG 2000 Image Compression (Lossless Only)	1.2.840.10008.1.2.4.90		
		JPEG 2000 Image Compression	1.2.840.10008.1.2.4.91		
		JPEG Baseline Lossy 8-Bit	1.2.840.10008.1.2.4.50		
		JPEG Baseline Lossy 12-Bit	1.2.840.10008.1.2.4.51		
		JPEG Lossless Hierarch., First-order prediction	1.2.840.10008.1.2.4.70		
Basic Text SR	1.2.840.10008.5.1.4.1.1.8 8.11	Explicit VR Little Endian	1.2.840.10008.1.2.1	SCP	None
		Implicit VR Little Endian	1.2.840.10008.1.2		
		Explicit VR Big Endian	1.2.840.10008.1.2.2		
Mammography CAD SR	1.2.840.10008.5.1.4.1.1.8 8.50	Explicit VR Little Endian	1.2.840.10008.1.2.1	SCP	None
		Implicit VR Little Endian	1.2.840.10008.1.2		
		Explicit VR Big Endian	1.2.840.10008.1.2.2		
Encapsulated PDF Storage	1.2.840.10008.5.1.4.1.1.104.1	Explicit VR Little Endian	1.2.840.10008.1.2.1	SCP	None
		Implicit VR Little Endian	1.2.840.10008.1.2		
		Explicit VR Big Endian	1.2.840.10008.1.2.2		
Storage Commitment Push Model	1.2.840.10008.1.20.1	Explicit VR Little Endian	1.2.840.10008.1.2.1	SCU	None
		Implicit VR Little Endian	1.2.840.10008.1.2		
		Explicit VR Big Endian	1.2.840.10008.1.2.2		

Note: The Storage Commitment Provider initiating the association must use the role selection negotiation.

2.3.3.1.2.1 SOP Specific DICOM Conformance Statement for Storage SOP Classes

Following are the status codes the application may send back to the SCU Equipment after performing the requested Storage:

Service Status	Status Codes	Further Meaning	Application Behaviour When receiving Status Codes	Related Fields Processed if received
----------------	--------------	-----------------	---	--------------------------------------



Refused	A7xx	Out of resources	Indicated that there was not enough space or some other internal resource to store the image.	(0000,0902)
Error	0110	Processing Failure	Indicates that an internal system call has failed while processing the image	(0000,0902)
Success	0000			None

Note: Images sent in JPEG-2000, JPEG 2000 (Lossless Only), JPEG Baseline Lossy 8-Bit, JPEG Baseline Lossy 12-Bit, JPEG Lossless Hierarch., First-order prediction transfer syntax are decompressed by the NETWORK AE by default.

Note: Correct transmission of images through the Network AE by itself does not guarantee an image to be visualized on MammoWorkstation. All modalities are managed by the database.

For images to be displayed on MammoWorkstation they must meet certain minimum DICOM requirements. Functionality, such as automatically sorted images in the Patient Folder can only be guaranteed when these DICOM requirements are met:

1. Accepted Modalities for X-ray mammography

General Series Module Attributes (C.7.3.1)

Tag	Attribute Name	Accepted Values
0008,0060	Modality	MG CR* OT
0018,0015	Body Part Examined	BREAST

Images from other modalities are not interpreted as X-Ray Mammography. Instead they appear in the “Additional Series” list.

*For Fuji-CR there is a “MG-Option” which is to be used instead of CR.

2. Sorted Patient Folder

General Image Module Attributes (C.7.6.1)

Tag	Attribute Name	Accepted Values
0020,0020	Patient Orientation	P,FL A,FR P,L A,R P,F A,F

Mammography Image Module Attributes (C.8.11.7)

Tag	Attribute Name	Accepted Values
-----	----------------	-----------------



0020,0062	Image Laterality	R L
-----------	------------------	-------

CR-Image (C.8.1.2), X-Ray Image (C.8.7.1), DX Positioning Module Attr. (C.8.11.5)

Tag	Attribute Name	Accepted Values
0018,5101	View Position	CC MLO ML

Remark: For Agfa CR also the Series Description (0008,103E), General Series Module Attributes C.7-5) is interpreted.

Translation table for standard orientations:

Patient Orientation	Image Laterality	View Position
P, L	R	CC
A, R	L	CC
P, FL	R	MLO
A, FR	L	MLO
P, F	R	ML
A, F	L	ML

MammoWorkstation supports additional values for patient orientation with automatic flip and rotation of images for proper display.

If Patient Orientation is missing, the „View Code Sequence“ is used according to SNM3:

Code Value	Code Meaning	equivalent View Position	View Position in Workstation
R-10242	cranio caudal	CC	CC
R-10224	medio-lateral	ML	ML
R-10226	medio-lateral oblique	MLO	MLO
R-10244	caudo-cranial	FB	FB
R-10228	latero-medial	LM	LM
R-10230	latero-medial oblique	LMO	LMO



R-102D0	superolateral to inferomedial oblique	SIO	SIO
R-102CF (deprecated)	cranio-caudal exaggerated	XCC	XCC
R-1024A	cranio-caudal exaggerated laterally	XCCL	XCCL
R-1024B	cranio-caudal exaggerated medially	XCCM	XCCM
Y-X1770 (retired)	cranio-caudal exaggerated laterally	XCCL	XCCL
Y-X1771 (retired)	cranio-caudal exaggerated medially	XCCM	XCCM

The following "View Modifier Sequence" values are supported:

Code Value	Code Meaning	Image type in Workstation (ACR-BIRADS code in parentheses)
R-102D6	Magnification	Magnification (M)
R-102D7	Spot Compression	Spot Compression (S)
R-102D2	Cleavage	Cleavage (CV)
R-102D5	Implant Displaced	Implant Displaced (ID)
R-102D1	Axillary tail	Axillary tail (AT)
R-102D3	Rolled lateral	Rolled lateral (RL)
R-102D4	Rolled medial	Rolled medial (RM)
R-102CA	Rolled inferior	Rolled inferior (RI)
R-102C9	Rolled superior	Rolled superior (RS)
R-102C2	Tangential view	Tangential (TAN)

If none of the above values are set, then the MammoWorkstation cannot sort the images according to laterality and view position. The images will then be displayed unsorted in the Patient Folder.

3. Measurements and Scales, Quadrant-Zoom, DICOM-Print

First, the software checks if Imager Pixel Spacing and "Estimated Radiographic Magnification Factor":



Tag	Attribute Name	Affected DICOM Modules
0018,1164	Imager Pixel Spacing	CR-Image (C.8.1.2) X-Ray Acquisition (C.8.7.2) DX Detector (C.8.11.4)

If the "Estimated Radiographic Magnification Factor" is present and the above Imager Pixel Spacing, then all measurement calculations are adjusted by the value of the magnification factor:

Tag	Attribute Name	Affected DICOM Modules
0008,1114	Estimated Radiographic Magnification Factor	DX Positioning (C.8.11.5)

Note: For modalities **GE Senographe 2000D, DS, Essential** the magnification factor is already pre-applied to the Imager Pixel Spacing. The MammoWorkstation will therefore ignore this factor on images of these modalities for measurements.

If either "Imager Pixel Spacing" or the "Estimated Radiographic Magnification Factor" is missing and the below value of "Pixel Spacing" is present, it will be used instead:

Tag	Attribute Name	Affected DICOM Modules
0028,0030	Pixel Spacing	Image Plane (C.7.6.2) CR-Image (C.8.1.2) DX Detector (C.8.11.4)

If both "Imager Pixel Spacing" and "Pixel Spacing" are missing (or the resulting value is zero), then the 1:1cm display function will be deactivated. Furthermore, the scales in the Magnifying Glass and measurement function will be deactivated. DICOM Print functionality will be limited as well, i.e. true-size printing is not available then.

4. Interpretation of different values for row Pixel Spacing and column Pixel Spacing

Although row and column values of Imager Pixel Spacing (or Pixel Spacing) might be different (i.e. anisotropic) the MammoWorkstation expects both row Pixel Spacing and column Pixel spacing to be identical. Anisotropic images are therefore neither supported for viewing, nor for any measurement, nor for any other activity. This implementation is in line with the DICOM Correction Item CP-626 "Clarify pixel spacing row and column order".

5. Windowing

VOI LUT Module Attributes (C.11.2)

0028,1050 Window Center

0028,1051 Window Width

0028,1055 Window Center & Width Explanation

0028,1056 VOI LUT Function



0028,3010 VOI LUT Sequence

If a VOI LUT Sequence is present, it overrides the Window Function attributes.

If none of these values are set, the MammoWorkstation applies default values.

In addition, the VOI LUT Module is supported, cf. Appendix

2.3.3.1.3 VOI LUT Attributes Module

Attribute name	Tag	Attribute Description
VOI LUT Sequence	(0028,3010)	Defines a sequence of VOI LUTs.
> LUT Descriptor	(0028,3002)	Specifies the format of the LUT Data in this Sequence. See C.11.2.1.1 for further explanation. Required if the VOI LUT Sequence (0028,3010) is sent..
> LUT Explanation	(0028,3003)	Free form explanation of the meaning of the LUT
> LUT Data	(0028,3006)	LUT Data in this Sequence. Required if the VOI LUT Sequence (0028,3010) is sent.
Window Center	(0028,1050)	Window Center for display. See C.11.2.1.2 for further explanation
Window Width	(0028,1051)	Window Width for display. See C.11.2.1.2 for further explanation. Required if Window Center (0028,1050) is sent.
Window Center & Width Explanation	(0028,1055)	Free form explanation of the meaning of the Window Center and Width. Multiple values correspond to multiple Window
VOI LUT Function	(0028,1056)	Defines a VOI LUT function to apply to the values of Window Center and Window Width. Supported values: LINEAR (default) SIGMOID

2.3.3.1.4 Presentation Context Acceptance Criterion

Only known SOP Classes are accepted.

2.3.3.1.5 Storage commitment (storage commitment broker)

Storage commit requests for DICOM data sent to the MammoWorkstation by a sending AE are never processed by the MammoWorkstation itself, because the MammoWorkstation is not a persistent DICOM archive. Instead, the MammoWorkstation forwards these storage commitment requests to another receiving AE, which is configured as DICOM partner in the MammoWorkstation. This functionality of the NETWORK AE is called "storage commitment broker". The declaration of remote hosts supporting storage commitment is done through a specific menu (DICOM Partner menu).



The storage commitment is initiated if the sending AE requests a storage commitment from the NETWORK AE for images sent to the MammoWorkstation. The MammoWorkstation forwards this storage commitment to the configured receiving SCP.

When a successful commitment for some received documents is received, the NETWORK AE forwards the result to the sending AE, which is the originator of the images and the originator storage commitment.

The NETWORK AE initiates the following operations:

- Initiate a DICOM association to ask a remote host (Storage Commitment SCP) on specific images.

The NETWORK AE waits for association requests from remote storage commitment AE:

- Answer to DICOM associations transmitting Storage Commitment Notification (N-EVENT-REPORT)

2.4 PRINT AE

The Print AE provides standard conformance to the following DICOM V3.0 Service Object Pair (SOP) Class as DICOM Print SCU.

SOP Class Name	SOP Class UID
Basic Grayscale Print Management (META)	1.2.840.10008.5.1.1.9
Basic Film Session	1.2.840.10008.5.1.1.1
Basic Film Box	1.2.840.10008.5.1.1.2
Basic Grayscale Image Box	1.2.840.10008.5.1.1.4
Printer	1.2.840.10008.5.1.1.16
Presentation LUT	1.2.840.10008.5.1.1.23

2.4.1 Association establishment policies

2.4.1.1 General

The PRINT AE will initiate an association as DICOM Basic Print Class User.

The PDU size is 64234 bytes.

The Print Management Service Class does not support extended negotiation.

2.4.1.2 Number of associations

The PRINT AE supports only one association at a time. Requests are internally queued.



2.4.1.3 Asynchronous nature

The PRINT AE does not support asynchronous communication (multiple outstanding transactions over a single association).

2.4.1.4 Implementation identifying information

The Implementation Class Unique Identifier (UID) for the PRINT AE is:

Implementation UID	1.2.276.0.69.25.1.1
--------------------	---------------------

2.4.2 Association initiation by real-world activity

There is only one real world activity that initiates association for the PRINT AE:

1. Manual Image Print

The PRINT AE initiates one association with the selected remote DICOM Printer. No other association can be opened by the PRINT AE while the current association is active.

2.4.2.1 Real-world activity: “Manual Image Print”

2.4.2.1.1 Associated real-world activity for “Manual Image Print”

The user has the possibility to select an image in the Viewer and to print it through the “Print” function of the Print Manager. The Print Manager allows to manipulate some print parameters like the number of copies. When the user requests to print by pushing the “Print” button, the PRINT AE tries to establish the association with the requested printer and sends the images to printing.

2.4.2.1.2 Proposed Presentation contexts for “Manual Image Print”

Presentation Context Table – Proposed					
Abstract Syntax		Transfer Syntax		Role	Extended Negotiation
Name	UID	Name List	UID List		
Basic Grayscale Print Management Meta SOP Class	1.2.840.10008.5.1.1.9	Implicit VR Little Endian	1.2.840.10008.1.2	SCU	None
Presentation LUT	1.2.840.10008.5.1.1.23	Implicit VR Little Endian	1.2.840.10008.1.2	SCU	None

2.4.2.1.2.1 SOP Specific DICOM Conformance Statement for Print Management SOP Class

For each of the supported Print Management SOP and Meta SOP Classes, the optional attributed and service elements supported, the valid range of values for mandatory and optional attributes, and the status code behavior are described in Annex.



2.4.3 Association Acceptance Policy

PRINT AE does not respond to attempts by a remote AE to open an association.



2.5 MEDIA INTERCHANGE AE

2.5.1 Introduction

The Media Interchange AE provides standard conformance to DICOM V3.0 Media Interchange. Media created complies with the IHE PDI profile. Service Object Pair (SOP) Class as File-Set Creator (FSC) and File-Set Reader (FSR).

MammoWorkstation provides capabilities to DICOM interchange on CD-Rs (Compact Disc-Recordable), and on CDROMs (Compact Disc Read Only Memory).

2.5.2 Implementation Model

2.5.2.1 Application Data Flow Diagram

Please confer to Figure 1.1.1.

2.5.3 Functional definition of the Media Interchange AE

The Media Interchange AE supports the following features:

- access to patient demographics and pixel data in the local database.
- generate a DICOM File Set (FSC) on the local disk to be burned on a CD-R or DVD-R.
- read a DICOM File Set (FSR) on a CD-R/DVD/CDROM or from local disk for showing the media directory.
- read a DICOM File Set (FSR) on a CD-R/ DVD/CDROM or from local disk for importing the selected studies.

2.5.4 File Meta Information For Implementation Class And Version

File Meta Information Version	1
Implementation UID	1.2.276.0.69.25.1.1

2.5.5 AE SPECIFICATION

2.5.5.1 Media Interchange AE Specification

The DICOM CDR/CDROM SERVER Application Entity provides standard conformance to

DICOM Interchange Option of the Media Storage Service Class. The application Profiles and roles are listed below.

Application Profiles Supported	Real World Activity	Role	SC Option
STD-AP	CD-R burning	FSC	Interchange



	CD-R/CDROM reading	FSR	Interchange
--	--------------------	-----	-------------

2.5.5.2 File Meta Information for the Application Entity

File Meta Information Version	1
-------------------------------	---

2.5.5.3 Real world activities for the Media Interchange AE

2.5.5.3.1 Display Directory

2.5.5.3.2 Real world activity “CD-R burning”

In this fashion the Media Interchange AE acts as File-Set Creator (FSC) using the Interchange Option. After selecting at least one study the user can create a new DICOM file set by first clicking the DICOM CD button afterwards the Burn CD button. The Directory Information Module is always created and the Basic Directory IOD Definition Model is made from standard tags either of type 1 or 2 and from the Study Instance UID tag which is of type 1C.

2.5.5.3.2.1 Options

Following are the SOP Classes supported by the “CD-burning” activity:

Abstract Syntax		Transfer Syntax		Role
Name	UID	Name List	UID List	
Computed Radiography Image Storage	1.2.840.10008.5.1.4.1.1.1	Explicit VR Little Endian	1.2.840.10008.1.2.1	FSC
Digital X-Ray Image Storage – For Presentation	1.2.840.10008.5.1.4.1.1.1.1			
Digital X-Ray Image Storage – For Processing	1.2.840.10008.5.1.4.1.1.1.1.1			
Digital Mammography X-Ray Image Storage – For Presentation	1.2.840.10008.5.1.4.1.1.1.2			
Digital Mammography X-Ray Image Storage – For Processing	1.2.840.10008.5.1.4.1.1.1.2.1			
CT Image Storage	1.2.840.10008.5.1.4.1.1.2			
MR Image Storage	1.2.840.10008.5.1.4.1.1.4			
US Image Storage	1.2.840.10008.5.1.4.1.1.6.1			



Secondary Capture Image Storage	1.2.840.10008.5.1.4.1.1.7			
Grayscale Softcopy Presentation State Storage	1.2.840.10008.5.1.4.1.1.11.1			
Basic Text SR	1.2.840.10008.5.1.4.1.1.88.11			
Mammography CAD SR	1.2.840.10008.5.1.4.1.1.88.50			
Encapsulated PDF Storage	1.2.840.10008.5.1.4.1.1.104.1			

2.5.5.3.3 Real world activity “CD-R/CDROM reading”

In this fashion the Media Interchange acts as File-Set Reader (FSR) using the Interchange Option and the application entity supports any Basic Directory IOD made from tags either of type 1 (1C) or 2.

The user can browse a DICOM file set from a CD-R/CDROM by first clicking the DICOM CD button and then the Import CD button. At this point a new Dialog Window is displayed on the screen where there is the list of the studies burned on the CD-R/CDROM: the user may select one or more studies to import into the local database by clicking the Import Studies button.

2.5.5.3.3.1 Options

Following are the SOP Classes supported by the “CD-R/CDROM reading” activity:

Abstract Syntax		Transfer Syntax		Role
Name	UID	Name List	UID List	
Computed Radiography Image Storage	1.2.840.10008.5.1.4.1.1.1	Explicit VR Little Endian	1.2.840.10008.1.2.1	FSR
Digital X-Ray Image Storage – For Presentation	1.2.840.10008.5.1.4.1.1.1.1			
Digital X-Ray Image Storage – For Processing	1.2.840.10008.5.1.4.1.1.1.1.1			
Digital Mammography X-Ray Image Storage – For Presentation	1.2.840.10008.5.1.4.1.1.1.2			
Digital Mammography X-Ray Image Storage – For Processing	1.2.840.10008.5.1.4.1.1.1.2.1			
CT Image Storage	1.2.840.10008.5.1.4.1.1.2			



MR Image Storage	1.2.840.10008.5.1.4.1.1.4			
US Image Storage	1.2.840.10008.5.1.4.1.1.6.1			
Secondary Capture Image Storage	1.2.840.10008.5.1.4.1.1.7			
Grayscale Softcopy Presentation State Storage	1.2.840.10008.5.1.4.1.1.11.1			
Basic Text SR	1.2.840.10008.5.1.4.1.1.88.11			
Mammography CAD SR	1.2.840.10008.5.1.4.1.1.88.50			
Encapsulated PDF Storage	1.2.840.10008.5.1.4.1.1.104.1			

2.5.6 Augmented And Private Application Profiles

No augmented or private application profiles are implemented by this application entity.

2.5.7 Extensions, Specializations And Privatizations Of SOP Classes And Transfer Syntaxes

No extensions, specialization or privatizations of SOP classes or transfer syntaxes are implemented by this application entity



3 Communication profiles

3.1 Supported Communication Stacks

MammoWorkstation provides DICOM V3.0 TCP/IP Network Communication Support as defined in PS 3.8.

3.2 TCP/IP Stack

MammoWorkstation uses the MergeCOM-3 Advanced DICOM Tool Kit to communicate over the TCP/IP protocol stack on any physical interconnection media supporting the TCP/IP stack. The tool kit inherits the TCP/IP stack from the operating system upon which it executes.

3.2.1 Physical Media Support

The MammoWorkstation is indifferent to the physical medium over which TCP/IP executes; it inherits this from the operating system on which it exists.

3.2.2 Internet Protocol Version 6 (IPv6)

Internet Protocol Version 6 (IPv6) is not supported by the MammoWorkstation.

4 Extensions/specializations/privatizations

4.1 Standard extended/specialized/private SOPs

None supported.

4.2 Private Transfer Syntaxes

None supported.

5 Configuration

The DICOM configuration is set by a MammoWorkstation configuration dialog available from the Login Screen. The operator may set the following parameters:

- MammoWorkstation AE Title and receive SCP Port
- Add or remove Sources (DICOM Query/Retrieve, DICOM Media)
- Add or remove Destinations (DICOM Store)
- Configure DICOM Modality Worklist SCU
- Configure remote DICOM Printer



Note: All configurations must be performed by a GE Field Engineer or a certified partner.

5.1 AE title/presentation address mapping

Before communicating with a remote AE (WORKLIST AE, NETWORK AE, PRINT AE) the operator must register it by using the MammoWorkstation DICOM configuration. This task requires to specify the following information:

- Remote AE Title
- Remote IP Address
- Remote TCP/IP Port Number

This information is used to communicate over the TCP/IP protocol stack.

6 Support of extended character sets

Not supported.



7 Annex: Structured Reporting SOP Class Definition

7.1 Introduction

This section of the DICOM Conformance Statement specifies the supported DICOM Structured Report SOP Classes, the optional attributes supported, and the valid range of values for mandatory and optional attributes.

7.2 Structured Report SOP Class Definitions

Image Diagnost follows DICOM Standard Part 16 und Supplement 79 for encoding for Mammography Structured Reports.

7.2.1 Templates

7.2.1.1 TID 4200 BREAST IMAGING REPORT

	NL	Rel with Parent	VT	Concept Name	VM	Req Type	Condition	Value Constraint	Set
1			CONTAINER	EV (111400, DCM, "Breast Imaging Report")	1	M			
2	>	HAS CONCEPT MOD	INCLUDE	DTID (1204) Language of Content Item and Descendants	1	M			
3	>	CONTAINS	INCLUDE	DTID (4202) Breast Imaging Report Narrative	1	M			
4	>	CONTAINS	INCLUDE	DTID (4208) Breast Imaging Report Supplementary Data	1	U			

7.2.1.2 TID 1204 LANGUAGE OF CONTENT ITEM AND DESCENDANTS

	NL	Rel with Parent	VT	Concept Name	VM	Req Type	Condition	Value Set Constraint
1		CONCEPT MOD	CODE	EV (121049,DCM,"Language of Content Item and Descendants"	1	M		DCID(5000) Image Diagnost wert ist: (RFC 3066, de-1996)
2		CONCEPT MOD	CODE	(121046,DCM,"Country of Language")	1	U		DCID(5001)



7.2.1.3 TID 4202 BREAST IMAGING REPORT NARRATIVE

	NL	Rel with Parent	VT	Concept Name	VM	Req Type	Condition	Value Constraint	Set
1			CONTAINER	EV (111412, DCM, "Narrative Summary")	1	M			
2	>	CONTAINS	CONTAINER	BCID (6052) Breast Imaging Report Section Title	1-n	M			
3	>>	HAS CONTEXT	OBS INCLUDE	DTID (1002) Observer Context	1-n	U	Not used by Image Diagnost		
4	>>	CONTAINS	TEXT	BCID (6053) Breast Imaging Report Elements	1	M			
5	>>>	INFERRED FROM	INCLUDE	DTID (350) References to Supporting Evidence	1	U	Not used by Image Diagnost		

7.2.1.4 TID 4208 BREAST IMAGING REPORT SUPPLEMENTARY DATA

	NL	Rel with Parent	VT	Concept Name	VM	Req Type	Condition	Value Set Constraint
1			CONTAINER EV	(111414, DCM, "Supplementary Data")	1	M		
2	>	CONTAINS	INCLUDE	DTID (4201) Breast Imaging Procedure Reported	1-n	M		
3	>	CONTAINS	CODE	EV (111403, DCM, "Baseline screening mammogram")	1	U	Not used by Image Diagnost	DCID (230) Yes-No
4	>	CONTAINS	CODE	EV (111404, DCM, "First mammogram ever")	1	U	Not used by Image Diagnost	DCID (230) Yes-No
5	>	CONTAINS	INCLUDE	DTID (4205) Breast Composition Section	1	U		
6	>	CONTAINS	INCLUDE	DTID (4206) Breast Imaging Report Finding Section	1-n	U		
7	>	CONTAINS	INCLUDE	DTID (4204) Breast Imaging Report Intervention Section	1-n	U	Not used by Image Diagnost	
8	>	CONTAINS	CONTAINER	EV (IDI, I-420010) Additional Information	1	U		Image Diagnost Extension
9	>>	CONTAINS	CODE	EV (IDI, I-604400) Image Quality Assessment	1-2	U		Image Diagnost Extension
10	>>	CONTAINS	CONTAINER EV	(111413, DCM, "Overall Assessment")	1	U		



11	>>>	CONTAINS	INCLUDE	DTID (4203) Breast Imaging Assessment	1	M		
12	>>	CONTAINS	CONTAINER	EV (IDI, I-420025) Conference Participants	1	U		Image Diagnost Extension
13	>>>	CONTAINS	CODE	EV (IDI, I-X990500) "First reading physician present"	1	U		Image Diagnost Extension
13	>>	CONTAINS	TEXT	(IDI, IMAGE_ITEM, "Image Item")	1-n	U		Image Diagnost Extension
14	>>	CONTAINS	TEXT	(121106, DCM, "Comment")	1	U		

7.2.1.5 TID 4205 BREAST COMPOSITION SECTION

	NL	Rel Parent	with	VT	Concept Name	VM	Req Type	Condition	Value Set Constraint
1				CONTAINER	EV (F-01710, SRT, "Breast composition")	1	M		
2	>	CONTAINS		INCLUDE	DTID (4201) Breast Imaging Procedure Reported	1-n	U	Not used by Image Diagnost	
3	>	CONTAINS		CODE	EV (F-01710, SRT, "Breast composition")	1-n	MC	At least one of row 3, 5 shall be present	DCID (6000) Overall Breast Composition
4	>>	HAS CONCEPT MOD		CODE	EV (G-C171, SRT, "Laterality") 1		M		DCID (6022) Side
5	>	CONTAINS		NUM	EV (111046, DCM, "Percent Glandular Tissue")	1-n	MC	At least one of row 3, 5 shall be present. Not used by Image Diagnost	UNITS = (% , UCUM "Percent")
6	>>	HAS CONCEPT MOD		CODE	EV (G-C171, SRT, "Laterality") 1		M	Not used by Image Diagnost	DCID (6022) Side

7.2.1.6 TID 4206 BREAST IMAGING REPORT FINDING SECTION

	NL	Rel Parent	with	VT	Concept Name	VM	Req Type	Condition	Value Set Constraint
1				CONTAINER	EV (121070, DCM, "Findings")	1	M		
2	>	HAS OBS CONTEXT		INCLUDE	DTID (1002) Observer Context	1-n	U	Not used by Image Diagnost	
3	>	CONTAINS		INCLUDE	DTID (4201) Breast Imaging Procedure Reported	1	M		



4	>	CONTAINS	CODE	EV (121071, DCM, "Finding")	1-n	M		DCID (6054) Breast Imaging Findings (includes CID 6016, 6057, 6064)
5	>>	HAS CONCEPT MOD	CODE	EV (111405, DCM, "Implant type") 1-n		UC	May be present if value of row 4 is (A-04010, SRT, "Implant") . Not used by Image Diagnost	DCID (6059) Breast Implant Types
6	>>	HAS PROPERTIES	INCLUDE	DTID (4203) Breast Imaging Assessment	1	U		
7	>>	HAS PROPERTIES	CODE	EV (111014, DCM, "Clockface or region")		U	only if „Region in one plane“ not set	DCID (6018) Clockface Location or Region
8	>>	HAS PROPERTIES	CODE	EV (111048, DCM, "Quadrant location")		U	Not used by Image Diagnost	DCID (6020) Quadrant Location
9	>>	HAS PROPERTIES	CODE	EV (I-601800, IDI, „Region in one plane“)	1	U	only if "Clockface or region" not set	Image Diagnost Extension
10	>>	HAS PROPERTIES	INCLUDE	DTID (1400) Linear Measurement	1-n	U		
11	>>	HAS PROPERTIES	INCLUDE	DTID (1401) Area Measurement	1-n	U	Not used by Image Diagnost	
12	>>	HAS PROPERTIES	INCLUDE	DTID (1402) Volume Measurement	1-n	U	Not used by Image Diagnost	
13	>>	HAS PROPERTIES	CODE	EV (111020, DCM, "Depth") 1		U	Not used by Image Diagnost	DCID (6024) Depth
14	>>	HAS PROPERTIES	CODE	EV (111035, DCM, "Lesion Density") 1		U	Not used by Image Diagnost	DCID (6008) Density Modifier
15	>>	HAS PROPERTIES	CODE	EV (M-020F9, SRT, "Shape") 1-n		U	Not used by Image Diagnost	DCID (6004) Mammography Characteristics of Shape
16	>>	HAS PROPERTIES	CODE	EV (111037, DCM, "Margins") 1-n		U	Not used by Image Diagnost	DCID (6006) Mammography Characteristics of Margin
17	>>	HAS PROPERTIES	CODE	EV (111009, DCM, "Calcification Type") 1-n		U	Not used by Image Diagnost	DCID (6010) Mammography Calcification Types
18	>>	HAS PROPERTIES	CODE	EV (111008, DCM, "Calcification Distribution") 1		U	Not used by Image Diagnost	DCID (6012) Calcification Distribution Modifier
19	>>	HAS PROPERTIES	NUM	EV (111038, DCM, "Number of calcifications") 1		U	Not used by Image Diagnost	UNITS = EV ({calcifications}, UCUM, "calcifications") Value = Integer 1 – n
20	>>	HAS PROPERTIES	CODE	EV (111407, DCM, "Implant finding") 1-n		U	Not used by Image Diagnost	DCID (6072) Breast Implant Findings



21	>>	HAS PROPERTIES	CODE	EV (G-C189, SRT, "Associated Finding") 1-n		U	Not used by Image Diagnost	DCID (6056) Associated Findings for Breast (includes CID 6015)
22	>>	HAS PROPERTIES	NUM	EV (111406, DCM, "Number of similar findings") 1		U	Not used by Image Diagnost	UNITS = EV ({findings}, UCUM, "findings") Value = Integer 2 - n
23	>>	HAS PROPERTIES	CODE	EV (F-01720, SRT, "Change since last mammogram") 1-n		U	Not used by Image Diagnost	DCID (6002) Change Since Last Mammogram or Prior Surgery
24	>>	INFERRED FROM	INCLUDE	DTID (350) References to Supporting Evidence	1	U	Not used by Image Diagnost	
25	>	CONTAINS	CODE	EV (IDI, I-420012) Finding Position in Report GUI	1-n	U		Image Diagnost Extension
26	>>	HAS PROPERTIES	TEXT	(IDI, MLO_POS_X, "MLO Position X")	1	U		Image Diagnost Extension
27	>	CONTAINS	INCLUDE	Apply DTID (4206) Breast Imaging Report Finding Section with modifications: 1. CONTAINER EV is (IDI, I-420011) Related Finding in Other Procedure 2. Laterality is not used	1-n	U		Image Diagnost Extension

7.2.1.7 TID 4201 BREAST IMAGING PROCEDURE REPORTED

	NL	Rel with Parent	VT	Concept Name	VM	Req Type	Condition	Value Set Constraint
1			CODE	EV (121058, DCM, "Procedure reported")	1	M		DCID (6050) Breast Procedure Reported Image Diagnost supported values: DCM 111408 Film Screen Mammography DCM 111409 Digital Mammography SRT P5-B8500 Ultrasonography of breast SRT P5-0900D MRI of breast
2	>	HAS CONCEPT MOD	CODE	EV (111464, DCM, "Procedure Modifier") 1-	n	U	Not used by Image Diagnost	DCID (6058) Procedure Modifiers for Breast



3	>	HAS CONCEPT MOD	CODE	EV (G-C171, SRT, "Laterality") 1		M		DCID (6022) Side
4	>	HAS PROPERTIES	CODE	EV (111401, DCM, "Reason for procedure") 1		U	Not used by Image Diagnost	DCID (6051) Breast Procedure Reason Image Diagnost supported values SRT R-42453 Screening DCM 111402 Clinical finding
5	>>	HAS CONCEPT MOD	CODE	EV (G-D709, SRT, "Relative time") 1		U	Not used by Image Diagnost	DCID (12102) Temporal Periods Relating To Procedure or Therapy
6	>>	HAS CONCEPT MOD	CODE	EV (111402, DCM, "Clinical Finding") 1-	n	UC	IFF row 4 value is "Clinical Finding". Not used by Image Diagnost	DCID (6055) Breast Clinical Finding or Indicated Problem
7	>>>	HAS PROPERTIES	CODE	EV (G-C171, SRT, "Laterality") 1		U	Not used by Image Diagnost	DCID (6022) Side
8	>	HAS PROPERTIES	DATE	EV (111060, DCM, "Study Date") 1		U	Not used by Image Diagnost	

7.2.1.8 TID 4203 BREAST IMAGING ASSESSMENT

NL	Rel with Parent	VT	Concept Name	VM	Req Type	Condition	Value Set Constraint
1		CODE	EV (111005, DCM, "Assessment Category")	1	M		DCID (6026) Mammography Assessment
2		CODE	EV (111053, DCM, "Recommended Follow-up")	1-n	U		BCID (6028) Mammography Recommended Follow-up
3	>	HAS CONCEPT MOD	CODE	EV (G-C171, SRT, "Laterality") 1		U	Not used by Image Diagnost DCID (6022) Side
4	>	HAS PROPERTIES	NUM	EV (111055, DCM, "Recommended Follow-up Interval") 1		U	Not used by Image Diagnost UNITS = DCID (6046) Units of Follow-up Interval Values = Integer ≥ 0, where 0 = immediate follow-up
5	>	HAS PROPERTIES	DATE	EV (111054, DCM, "Recommended Follow-up Date") 1		U	Not used by Image Diagnost
6	>	HAS PROPERTIES	INCLUDE	DTID (4207) Breast Imaging Pathology Results	1-n	U	Not used by Image Diagnost



7.2.1.9 TID 1400 LINEAR MEASUREMENT

	NL	Rel Parent	with	VT	Concept Name	VM	Req Type	Condition	Value Set Constraint
1				NUM	DCID (7470) "Linear Measurements" Image Diagnost is using SNM3 M-02550 Diameter	1	M		UNITS = DCID(7460) "Units of Linear Measurement"
2	>	INFERRED		SCoord	EV (121055,DCM, "Path")	1	UC	XOR Row 5. Not used by Image Diagnost	GRAPHIC TYPE = {POLYLINE, CIRCLE, ELLIPSE}
3	>>	RSELECTED		IMAGE		1	MC	XOR Row 4. Not used by Image Diagnost	
4	>>	SELECTED		IMAGE		1	MC	XOR Row 3. Not used by Image Diagnost	
5	>>	INFERRED		SCoord	SCoord EV (121230, DCM, "Path Vertex")	2-n	UC	XOR Row 2. Not used by Image Diagnost	GRAPHIC TYPE = {POINT}
6	>>	RSELECTED		IMAGE		1	MC	XOR Row 7. Not used by Image Diagnost	
7	>>	SELECTED		IMAGE		1	MC	XOR Row 6. Not used by Image Diagnost	



8 Annex: Print Management SOP Class Definition

8.1 Introduction

This section of the DICOM Conformance Statement specifies the supported Print Management SOP and Meta SOP Classes, the optional attributes and service elements supported, the valid range of values for mandatory and optional attributes, and the status code behavior.

8.2 Print Management SOP Class Definitions

8.2.1 Basic Film Session SOP Class

The DICOM Print SCU AE supports the N-CREATE DIMSE Service Element for the Basic Film Session SOP Class.

- The N-CREATE DIMSE Service element sent by the DICOM Print SCU AE requests the Remote DICOM Print SCP to create an instance of Basic Film Session.

8.2.1.1 IOD Description

8.2.1.1.1 IOD Modules

Module	Reference	Module Description
SOP Common	C.12.1	Contains SOP Common information
Basic Film Session Presentation Module	C.13.1	Contains Film Session presentations information
Basic Film Session Relationship	C.13.2	References to related SOPs

8.2.1.1.2 Basic Film Session Presentation Module

Attribute name	Tag	Attribute Description
Number of Copies	(2000,0010)	1 to 99; default is set to 1.
Print Priority	(2000,0020)	Depending on configuration of associated Remote DICOM printer. Supported values HIGH MED LOW
Film Session Label	(2000,0050)	Empty



Medium Type	(2000,0030)	Depending on configuration of associated Remote DICOM printer. Supported values PAPER CLEAR FILM BLUE FILM MAMMO CLEAR FILM MAMMO BLUE FILM
Film Destination	(2000,0040)	Empty

8.2.1.1.3 Basic Film Session Relationship Module

Attribute name	Tag	Attribute Description
Referenced Film Box Sequence	(2000,0500)	Empty
>Referenced SOP Class UID	(0008,1150)	Empty
>Referenced SOP Instance UID	(0008,1155)	Empty

8.2.1.2 DIMSE Service Group

DIMSE Service Element	Usage SCU
N-CREATE	M
N-SET	Not used
N-DELETE	Not used
N-ACTION	Not used

8.2.1.2.1 N-CREATE

8.2.1.2.1.1 Attributes

Attribute name	Tag	Usage SCU
Number of Copies	(2000,0010)	Used



Print Priority	(2000,0020)	Used
Medium Type	(2000,0030)	Used
Film Session Label	(2000,0050)	Not used
Memory Allocation	(2000,0060)	Not used
Film Destination	(2000,0040)	Not used

8.2.1.2.1.2 Status

Service Status	Status Codes	Further Meaning	Application Behaviour When receiving Status Codes
Warning	B600	Memory allocation not supported	Association is aborted
Success	0000	Film session successfully created	Next step describe in the sequencing of Real-World Activities paragraph is performed

8.2.1.2.1.3 Behavior

No specific behavior.

8.2.2 Basic Film Box SOP Class

The DICOM Print SCU AE supports the following DIMSE Service Element for the Basic Film Box SOP Class.

- The N-CREATE DIMSE Service element sent by the DICOM Print SCU AE requests the Remote DICOM Print SCP to create an instance of Basic Film Box
- The N-ACTION DIMSE Service element sent by the DICOM Print SCU AE requests the Remote DICOM Print SCP to print the Basic Film Box onto the hard copy printer.
- The N-DELETE DIMSE Service element sent by the DICOM Print SCU AE requests the Remote DICOM Print SCP to release the Basic Film Box instance.

8.2.2.1 IOD Description

8.2.2.1.1 IOD Modules

Module	Reference	Module Description
SOP Common	C.12.1	Contains SOP Common information
Basic Film Box Presentation Module	C.13.3	Contains Film Box presentations information
Basic Film Box Relationship	C.13.4	References to related SOPs



8.2.2.1.2 Basic Film Box Presentation Module

Attribute name	Tag	Attribute Description
Image Display Format	(2010,0010)	STANDARD\1,1 STANDARD\1,2 STANDARD\2,2 STANDARD\2,3 STANDARD\3,3 STANDARD\3,4 STANDARD\3,5 STANDARD\4,4 STANDARD\4,5 STANDARD\4,6
Annotation Display Format ID	(2010,0030)	Not sent.
Film Orientation	(2010,0040)	PORTRAIT LANDSCAPE



Film Size ID	(2010,0050)	8INX10IN 10INX12IN 10INX14IN 11INX14IN 14INX14IN 14INX17IN 24CMX24CM 24CMX30CM DIN A4 DIN A3
Border density	(2010,0100)	Not sent.
Empty Image Density	(2010,0110)	Not sent.
Magnification Type	(2010,0060)	Optional (value depending on configuration of associated Remote DICOM printer)
Smoothing Type	(2010,0080)	Optional (value depending on configuration of associated Remote DICOM printer)
Min Density	(2010,0120)	Optional (value depending on configuration of associated Remote DICOM printer)
Max Density	(2010,0130)	Optional (value depending on configuration of associated Remote DICOM printer)
Trim	(2010,0140)	Not sent.
Configuration Information	(2010,0150)	Optional (value depending on configuration of associated Remote DICOM printer)
Requested Resolution ID	(2020,0050)	Not sent.

8.2.2.1.3 Basic Film Box Relationship Module

Attribute name	Tag	Attribute Description
----------------	-----	-----------------------



Referenced Film Session Sequence	(2010,0500)	Empty
>Referenced SOP Class UID	(0008,1150)	Empty
>Referenced SOP Instance UID	(0008,1155)	Empty
Referenced Image Box Sequence	(2010,0510)	Empty
>Referenced SOP Class UID	(0008,1150)	Empty
>Referenced SOP Instance UID	(0008,1155)	Empty
Referenced Basic Annotation Box Sequence	(2010,0520)	Empty
>Referenced SOP Class UID	(0008,1150)	Empty
>Referenced SOP Instance UID	(0008,1155)	Empty
Referenced Presentation LUT Sequence	(2050,0500)	Used if Presentation LUT SOP was created
>Referenced SOP Class UID	(0008,1150)	Used if Presentation LUT SOP was created
>Referenced SOP Instance UID	(0008,1155)	Used if Presentation LUT SOP was created

8.2.2.2 DIMSE Service Group

DIMSE Service Element	Usage SCU
N-CREATE	M
N-ACTION	M
N-DELETE	Used

8.2.2.2.1 N-CREATE

8.2.2.2.1.1 Attributes

Attribute name	Tag	Usage SCU
Image Display Format	(2010,0010)	M
Referenced Film Session Sequence	(2010,0500)	M



>Referenced SOP Class UID	(0008,1150)	M
>Referenced SOP Instance UID	(0008,1155)	M
Referenced Image Box Sequence	(2010,0510)	Not used
>Referenced SOP Class UID	(0008,1150)	Not used
>Referenced SOP Instance UID	(0008,1155)	Not used
Referenced Basic Annotation Box Sequence	(2010,0520)	Not used
>Referenced SOP Class UID	(0008,1150)	Not used
>Referenced SOP Instance UID	(0008,1155)	Not used
Film Orientation	(2010,0040)	Not used
Film Size ID	(2010,0050)	Used
Magnification Type	(2010,0060)	Used, not sent if empty
Smoothing Type	(2010,0080)	Used, not sent if empty
Min Density	(2010,0120)	Used, not sent if empty
Max Density	(2010,0130)	Used, not sent if empty
Configuration Information	(2010,0150)	Used, not sent if empty
Referenced Presentation LUT Sequence	(2050,0500)	Used
>Referenced SOP Class UID	(0008,1150)	Used
>Referenced SOP Instance UID	(0008,1155)	Used
Annotation Display Format ID	(2010,0030)	Not used
Smoothing Type	(2010,0080)	Not used
Border Density	(2010,0100)	Not used
Empty Image Density	(2010,0110)	Not used
Min Density	(2010,0120)	Not used
Trim	(2010,0140)	Not used



Requested Resolution ID	(2020,0050)	Not used
-------------------------	-------------	----------

8.2.2.2.1.2 Status

There are no specific status codes.

8.2.2.2.1.3 Behavior

There is no specific behavior.

8.2.2.2.2 N-DELETE

8.2.2.2.2.1 Behavior

The SCU uses the N-DELETE to request the SCP to delete the Basic Film Box SOP Instance hierarchy.

8.2.2.2.3 N-ACTION

N-ACTION is used to print the current film of the film session.

8.2.2.2.3.1 Attributes

Action Type Name	Action Type ID	Attribute	Tag	Usage SCU
Print	1	Referenced Print Job Sequence	(2100,0500)	Not used
		>Referenced SOP Class UID	(0008,1150)	Not used
		>Referenced SOP Instance UID	(0008,1155)	Not used

8.2.2.2.3.2 Status

Service Status	Status Codes	Further Meaning	Application Behaviour When receiving Status Codes
Success	0000	Film accepted for printing.	Next Step described in the sequencing of Real-World Activities paragraph is performed
Warning	B603	Film Box SOP Instance hierarchy does not contain Image Box SOP Instances (empty page)	This case should not happen. This warning is considered as an error. Association is aborted.
Failure	C602	Unable to create Print Job SOP Instance; print queue is full	Appropriate message is returned to the user. Association is aborted.



	C604	Image position collision : multiple images assigned to single image position	Appropriate message is returned to the user. Association is aborted.
	C603	Image size is larger than image box size (by using the specified magnification value)	Appropriate message is returned to the user. Association is aborted.

8.2.2.2.3.3 Behavior

SCU uses the N-ACTION to request the SCP to print one or more copies of a single film of the film session.

8.2.3 Image Box SOP Class

8.2.3.1 Basic Grayscale Image Box SOP Class

The DICOM Print SCU AE supports the following DIMSE Service Element for the Basic Grayscale Image Box SOP Class.

- The N-SET DIMSE Service element sent by the DICOM Print SCU AE requests the Remote DICOM Print SCP to set the attributes of the Basic Grayscale Image Box Instance.

8.2.3.1.1 IOD Description

8.2.3.1.1.1 IOD Modules

Module	Reference	Module Description
SOP Common	C.12.1	Contains SOP Common information
Image Box Presentation Module	C.13.5	Contains Image Box presentations information

8.2.3.1.1.2 Image Box Pixel Presentation Module

Attribute name	Tag	Attribute Description
Image Position	(2020,0010)	Position ID for composite printing (1 .. 4). Default is 1 for single-image printing.
Polarity	(2020,0020)	NORMAL = pixels shall be printed as specified by the Photometric Interpretation (0028,0004) REVERSE = pixels shall be printed with the opposite polarity as specified by the Photometric Interpretation (0028,0004)
Magnification Type	(2010,0060)	Not Sent.



Smoothing Type	(2010,0080)	Not Sent
Configuration Information	(2010,0150)	Not Sent
Requested Image Size	(2020,0030)	Set to Original Film Size if Requested Resolution ID is not set.
Preformatted Grayscale Image Sequence	(2020,0110)	This sequence is always included if the Image Box is a Basic Grayscale Image Box
>Samples Per Pixel	(0028,0002)	1
>Photometric Interpretation	(0028,0004)	Set from original image. Supported values MONOCHROME1 MONOCHROME2
>Rows	(0028,0010)	Height of preformatted pixmap in pixel
>Columns	(0028,0011)	Width of preformatted pixmap in pixel
>Pixel Aspect Ratio	(0028,0034)	1\1
>Bits Allocated	(0028,0100)	8 or 16, depending on source image data
>Bits Stored	(0028,0101)	8 or 12, depending on source image data. If depth of source image differs from 8 or 12 the image data will be re-quantized accordingly
>High Bit	(0028,0102)	Set to Bits Stored -1 (i.e. 7 or 11)
>Pixel Representation	(0028,0103)	0 (Unsigned Integer)
>Pixel Data	(7FE0,0010)	

8.2.3.1.2 DIMSE Service Group

DIMSE Service Element	Usage SCU
N-SET	M

8.2.3.1.2.1 N-SET

8.2.3.1.2.1.1 Attributes

Attribute name	Tag	Usage SCU
----------------	-----	-----------



Image Position	(2020,0010)	M
Preformatted Grayscale Image Sequence	(2020,0110)	M
>Samples Per Pixel	(0028,0002)	M
>Photometric Interpretation	(0028,0004)	M
>Rows	(0028,0010)	M
>Columns	(0028,0011)	M
>Pixel Aspect Ratio	(0028,0034)	1/1
>Bits Allocated	(0028,0100)	M
>Bits Stored	(0028,0101)	M
>High Bit	(0028,0102)	M
>Pixel Representation	(0028,0103)	M
>Pixel Data	(7FE0,0010)	M
Polarity	(2020,0020)	Used
Referenced Overlay Sequence	(0008,1130)	Not used
>SOP Class UID	(0008,1150)	Not used
>SOP Instance UID	(0008,1155)	Not used
Configuration Information	(2010,0150)	Not used
Magnification Type	(2010,0060)	Not used
Smoothing Type	(2010,0080)	Not used
Requested Image Size	(2020,0030)	Used, not sent if empty

8.2.3.1.2.1.2 Status

Service Status	Status Codes	Further Meaning	Application Behavior When receiving Status Codes
----------------	--------------	-----------------	--



Failure	C605	Insufficient memory in printer to store the image.	Appropriate message is returned to the user. Association is aborted.
---------	------	--	---

8.2.3.1.2.1.3 Behavior

There is no specific behavior.

The SCU does not instruct the SCP to erase the image in the image position by setting a zero length and no value in the Attribute Pre-formatted Grayscale Image Sequence (2020,0110) or Pre-formatted Grayscale Image Sequence (2020,0111).

8.2.4 Printer SOP Class

The DICOM Print SCU AE supports the following DIMSE Service Element for the Basic Printer SOP Class.

- The N-EVENT_REPORT DIMSE Service element sent by the DICOM Print SCP to the local DICOM Print SCU AE. The DICOM Print SCU handles the Printer Status and Printer Status Info fields. All other received data are ignored.
- The N-GET DIMSE Service element sent by the DICOM Print SCU AE requests the Remote DICOM Print SCP to give information on the Remote DICOM Printer.

8.2.4.1 IOD Description

8.2.4.1.1 IOD Modules

Module	Reference	Module Description
SOP Common	C.12.1	Contains SOP Common information
Printer Module	C.13.9	Contains status information to monitor printer

8.2.4.1.2 Printer Module

Attribute name	Tag	Attribute Description
Printer Status	(2110,0010)	The behavior defined for the following term NORMAL : Association is established. FAILURE : Association is released. WARNING : Association is released.
Printer Status Info	(2110,0020)	Printer shall return value



Printer Name	(2110,0030)	Printer shall return value
Manufacturer	(0008,0070)	Printer shall return value
Manufacturer Model Name	(0008,1090)	Printer shall return value
Device Serial Number	(0018,1000)	Printer shall return value
Software Versions	(0018,1020)	Printer shall return value
Date Of Last Calibration	(0018,1200)	Printer shall return value
Time Of Last Calibration	(0018,1201)	Printer shall return value

8.2.4.2 DIMSE Service Group

DIMSE Service Element	Usage SCU
N-EVENT-REPORT	M
N-GET	U

8.2.4.2.1 N-EVENT-REPORT

8.2.4.2.1.1 Attributes

Event Type Name	Event Type ID	Attribute	Tag	Usage SCU
Normal	1			
Warning	2	Printer Name	(2110,0030)	
		Printer Status Info	(2110,0020)	
Failure	3	Printer Name	(2110,0030)	
		Printer Status Info	(2110,0020)	

8.2.4.2.1.2 Behavior

If Printer Status is FAILURE or WARNING

Signal print failure to the user

Else

Signal print success to the user



8.2.4.2.2 N-GET

8.2.4.2.2.1 Attributes

Attribute name	Tag	Attribute Description
Printer Status	(2110,0010)	Used
Printer Status Info	(2110,0020)	Used
Printer Name	(2110,0030)	Used
Manufacturer	(0008,0070)	Used
Manufacturer Model Name	(0008,1090)	Used
Device Serial Number	(0018,1000)	Used
Software Versions	(0018,1020)	Used
Date Of Last Calibration	(0018,1200)	Used
Time Of Last Calibration	(0018,1201)	Used

8.2.4.2.2.2 Behavior

If Printer Status is FAILURE

Signal print failure to the user

Else If Printer Status is WARNING and Printer Status Info is not equal to SUPPLY LOW or SUPPLY EMPTY

Signal print warning to the user

Else

Signal print success to the user

8.2.5 Presentation LUT SOP Class

The DICOM Print SCU AE supports the following DIMSE Service Element for the Presentation LUT SOP Class.

- The N-CREATE DIMSE Service element sent by the DICOM Print SCU AE requests the Remote DICOM Print SCP to create an Instance of Presentation LUT.



8.2.5.1 IOD Description

8.2.5.1.1 IOD Modules

Module	Reference	Module Description
SOP Common	C.12.1	Contains SOP Common information
Presentation LUT Module	C.11.4	Contains Image presentation LUT information

8.2.5.1.2 Presentation LUT Module

Attribute name	Tag	Attribute Description
Presentation LUT Sequence	(2050,0010)	Defines a sequence of Presentation LUTs.
>LUT Descriptor	(0028,3002)	Depends on original image pixel depth. For example a 12 bit image leads to 4096\0\12
>LUT Explanation	(0028,3003)	Set to "Printer_PLUT"
>LUT Data	(0028,3006)	LUT Data in this Sequence.
Presentation LUT Shape	(2050,0020)	Not sent.

8.2.5.2 DIMSE Service Group

DIMSE Service Element	Usage SCU
N-CREATE	M
N-DELETE	Not used

8.2.5.2.1 N-CREATE

8.2.5.2.1.1 Attributes

Attribute name	Tag	Usage SCU
Presentation LUT Sequence	(2050,0010)	Used
>LUT Descriptor	(0028,3002)	Used
>LUT Explanation	(0028,3003)	Used
>LUT Data	(0028,3006)	Used



Presentation LUT Shape	(2050,0020)	Not sent.
------------------------	-------------	-----------

8.2.5.2.1.2 Status

Service Status	Status Codes	Further Meaning	Application Behavior When receiving Status Codes
Warning	B605	Requested Min Density or Max Density outside of printer's operating range. The printer will use its respective minimum or maximum density value instead.	
Success	0000	Presentation LUT successfully created	

8.2.5.2.1.3 Behavior

There is no specific behavior.



9 Annex: JPEG Decompression

The MammoWorkstation supports the JPEG_BASELINE, JPEG_EXTENDED_2_4, JPEG_LOSSLESS_HIER_14, JPEG_2000, JPEG_2000_LOSSLESS_ONLY transfer syntaxes with photometric interpretations MONOCHROME1, MONOCHROME2, RGB, and YBR. There are limits on the performance. For JPEG_BASELINE, JPEG_EXTENDED_2_4, and JPEG_LOSSLESS_HIER_14, images can be compressed or decompressed at a maximum rate of 3 images (or frames) per second. For JPEG_2000 and JPEG_2000_LOSSLESS_ONLY, a dialog will be displayed on Windows each time the compressor or decompressor is used.

9.1 Transfer syntax UID 1.2.840.10008.1.2.4.50, JPEG Baseline (Process 1)

Default Transfer Syntax for Lossy JPEG 8 Bit Image Compression. Table 19 details the photometric interpretation and bit depths supported for this transfer syntax.

JPEG Baseline				
Photometric Interpretation	MONOCHROME1 MONOCHROME2	RGB	YBR_FULL_422 YBR_FULL	
Bits Stored	8	8	8	
Bits Allocated	8	8	8	
Samples Per Pixel	1	3	3	

Table 19: JPEG Baseline Supported Photometric Interpretations and Bit Depths

9.2 Transfer syntax UID 1.2.840.10008.1.2.4.51, JPEG Extended (Process 2 & 4)

Table 20 details the photometric interpretation and bit depths supported by the decompressor for this transfer syntax.

JPEG Extended (Process 2 & 4)					
Photometric Interpretation	MONOCHROME1 MONOCHROME2	RGB	YBR_FULL_422 YBR_FULL		
Bits Stored	8	10	12	8	8
Bits Allocated	8	16	16	8	8
Samples Per Pixel	1	1	1	3	2

Table 20: JPEG Extended Supported Photometric Interpretations and Bit Depths



9.3 Transfer syntax UID 1.2.840.10008.1.2.4.70, JPEG Lossless, Non-Hierarchical, First-Order Prediction (Process 14 [Selection Value 1])

Table 21 details the photometric interpretation and bit depths supported by the decompressor for this transfer syntax.

JPEG Lossless Non-Hierarchical Process 14			
Photometric Interpretation	MONOCHROME1 MONOCHROME2		RGB
Bits Stored	2 - 8	9 - 16	8
Bits Allocated	8	16	8
Samples Per Pixel	1	1	3

Table 21: JPEG Lossless Supported Photometric Interpretations and Bit Depths

New compressor decompressor allow lossless Jpeg compression/decompression of images with the number of bits stored between 2 and 16. Previously only images with bits stored of 8, 10, 12 and 16 were accepted.

9.4 Transfer Syntax UID 1.2.840.10008.1.2.4.91, JPEG 2000 Image Compression

Table 22 details the photometric interpretation and bit depths supported by the compressor and decompressor for this transfer syntax.

JPEG 2000 (When used for Lossy)						
Photometric Interpretation	MONOCHROME1 MONOCHROME2				YBR_ICT	RGB
Bits Stored	8	10	12	16	8	8
Bits Allocated	8	16	16	16	8	8
Samples per Pixel	1	1	1	1	3	3

Table 22: JPEG 2000 Lossy Supported Photometric Interpretations and Bit Depths

9.5 Transfer syntax UID 1.2.840.10008.1.2.4.90, JPEG 2000 Image Compression (Lossless Only)

Table 23 details the photometric interpretation and bit depths supported by the compressor and decompressor for this transfer syntax.



JPEG 2000 Lossless						
Photometric Interpretation	MONOCHROME1 MONOCHROME2				YBR_RCT	RGB
Bits Stored	8	10	12	16	8	8
Bits Allocated	8	16	16	16	8	8
Samples Per Pixel	1	1	1	1	3	3

Table 23: JPEG 2000 Lossless Supported Photometric Interpretations and Bit Depths

All data needs to be right justified, i.e. bit 0 contains data, but the highest bits may not. RGB and YBR must be non-planar (R1G1B1, R2G2B2, ... or Y1Y2B1R1, Y3Y4B3R3,...) JPEG_2000 / JPEG_2000_LOSSLESS_ONLY will cause a irreversible, or reversible color transformation when compressing RGB data. The Photometric Interpretation is changed from RGB to:

- YBR_ICT for lossy color transform for lossy compression
- YBR_RCT for lossless color transform for lossless compression

Similarly, on the decompression end, the Photometric Interpretation is changed back to RGB.